

Review Manager 5 Tutorial

for Diagnostic Test Accuracy reviews

Welcome to the RevMan Diagnostic Test Accuracy (DTA) Tutorial. This tutorial is designed to give Cochrane review authors an introduction to the process of writing a Cochrane systematic review of the accuracy of diagnostic tests using RevMan. The approximate time needed to complete this tutorial will be about two to three hours.

The tutorial does not provide all the information necessary to write a Cochrane DTA review, and does not address other types of Cochrane systematic reviews such as reviews of intervention studies or overviews of reviews. For detailed instructions on how to conduct a Cochrane systematic review of an intervention, please refer to the *Cochrane Handbook for Systematic Reviews of Interventions*, available from the RevMan Help menu. For detailed instructions on how to conduct a Cochrane DTA review, please refer to the *Cochrane Handbook for Systematic Reviews of Diagnostic Test Accuracy*, available (in due time) from the following website
<http://srtda.cochrane.org>.

Later parts of this tutorial rely on activities completed in the earlier sections, so we recommend you start working from the beginning. If you would like to skip ahead, click on the section you would like to read in the Contents list.

If you have technical questions regarding RevMan software, please contact the Cochrane Information Management System team at www.cc-ims.net/RevMan5. Other questions related to preparing protocols or reviews in RevMan should be directed to the relevant Cochrane Review Group. Please send your comments about the Tutorial itself to cesu@amc.uva.nl.

Note: The results reported in this tutorial are fictional and were created by the authors.

PART 1 – Getting started	3
Opening your review.....	3
Setting up user details	3
The RevMan 5 interface.....	4
Saving a review	5
Adding authors to a review	5
PART 2 – Working with text.....	7
Adding text to a review.....	7
Formatting text.....	7
Making changes to headings	7
Adding notes	8
Tracking changes in the text	9
Checking your spelling	10
PART 3 – Adding studies and references.....	12
Adding studies to a review.....	12
Adding references to a study	12
Editing an existing reference.....	14
Importing references from a text file.....	15
Moving references.....	15
Adding links to references in the text.....	16
Removing a link	17
Adding multiple links automatically	17
PART 4 - Tables.....	18
Characteristics of included studies table	18
Assessment of methodological quality table.....	18
Adding self-defined quality items	19
Characteristics of excluded studies	20
Additional tables	21
Adding Footnotes to tables.....	21
PART 5 - Data and analyses	23
Adding tests.....	23
Entering data by test.....	23
Entering data using the calculator	25
Entering data by study.....	26
Assessment of methodological quality for two tests	27
Covariates.....	28
Analyses	30
a. Changing Forest Plot properties	31
b. Changing SROC Graph properties	32
c. Entering Externally Calculated Parameters.....	33
Test Comparisons.....	33
Paired Test Comparisons.....	35
Investigating heterogeneity	35
Adding a figure to your review	37
Methodological quality graphs and summaries.....	37
PART 6 – Finishing the review	39
Summary of results tables	39
Adding an appendix	39
Validating a review	39
Submitting a review for editorial approval	40

PART 1 – Getting started

Opening your review

All Cochrane reviews are stored in Archie, The Cochrane Collaboration's online database. When you have registered a review with a Cochrane Review Group, they will set up a user account for you in Archie. To work on your review, you will need to find your review online in Archie and check it out into RevMan using the **Check Out**  button on the toolbar.


For this tutorial we will use a sample review called 'Laboratory tests for caffeine addiction in PhD students'. Instead of checking the review out of Archie, it has been provided as an additional file with RevMan 5.

1. From the **File** menu, select **Open**, or click the **Open** button  on the toolbar.
2. Find the folder where RevMan was installed on your computer. If you are using Windows and accepted the default installation location, this will be **C:\Program Files\Review Manager 5**.
3. Open the **doc** folder, and then the **tutorial** folder.
4. Select the file 'Laboratory tests for caffeine addiction in PhD students.rm5'.
5. Click **Open**.

Note: It is possible to create a new review in RevMan by selecting **New** from the **File** menu, or clicking the **New** button  on the toolbar. Be aware that if you create a review yourself, you cannot check it into Archie or use it for publication in *The Cochrane Library*.

Setting up user details


If you will be using RevMan regularly on the same computer, you may wish to enter your user details (such as your name). RevMan will attach these user details to any editing done using track changes so that co-authors and editors can easily identify who has made the changes. Open your user details by selecting **Preferences** from the **Tools** menu. Enter your name directly or click **Choose** to find your details in Archie.

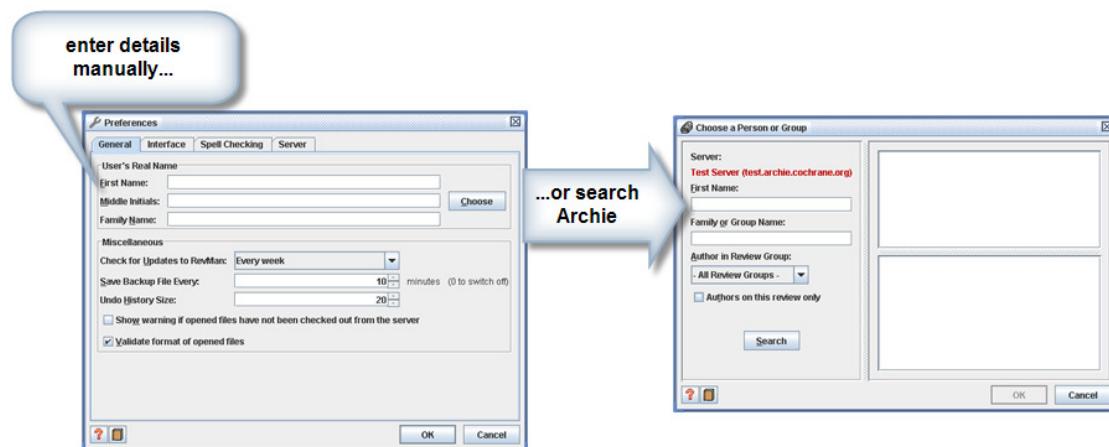

Note: You will need to set up RevMan access to Archie via the internet. In the **Preferences** window, open the **Connection** tab. Click the **Wizard** button to set up your connection, or enter your Archie user name and password directly and click the **Test** button to check the connection. You may need to enter proxy server details. Consult the RevMan 5 Help or your local IT support for more information.

The RevMan 5 interface

When a review is open in RevMan 5, you will see there are two separate window panes that interact with each other as you navigate through the different sections of a review. The pane on the left hand side shows the outline of a review. This is referred to as the *outline pane*. The area on the right hand side shows all the information contained in a review. This is referred to as the *content pane*. The *content pane* will usually show the text of your review, but can also contain additional tabs for results, figures, etc.

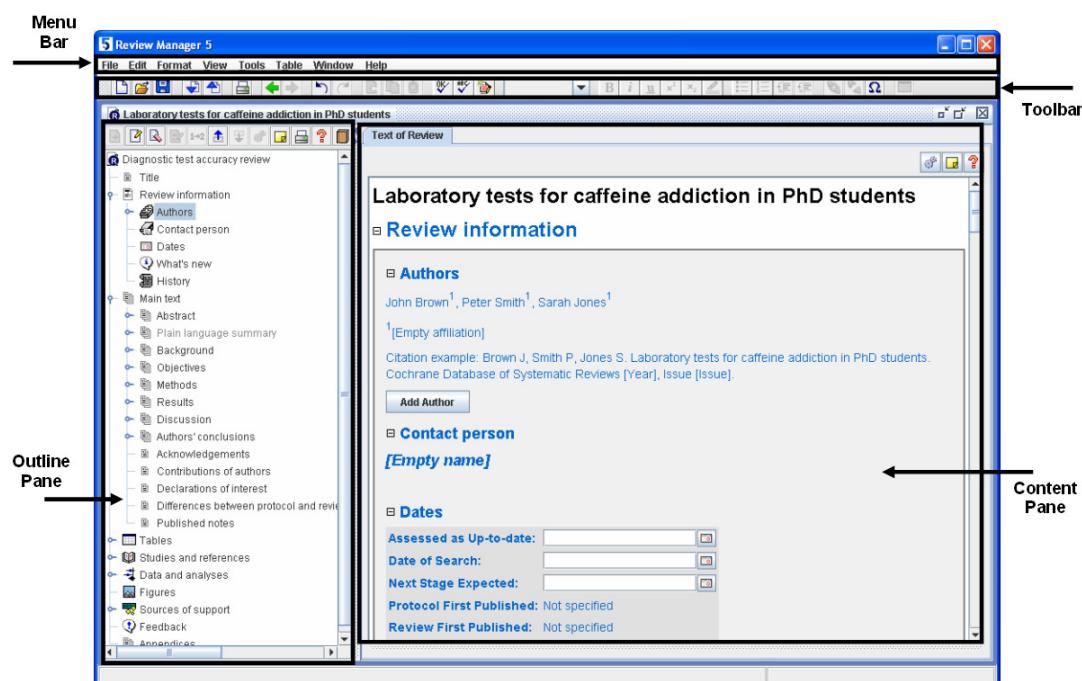

In the *outline pane*, you can use the key icon next to each heading to expand and collapse subheadings. In the *content pane*, use the symbol.

RevMan allows you to change the appearance of the interface by selecting **Preferences** from the **Tools** menu and going to the **Interface** tab. The default **Look and Feel** of the interface is **Metal**. All screen captures presented in this tutorial display the **Metal** look and feel. It is also possible to set language and spelling options in the **Preferences** window.

By default, your new review is at the **protocol stage**, indicating that you are planning the methods for a review. Sections of a review that are not required at the protocol stage (but will be needed later for your completed review) appear as inactive (light grey) in both the *outline pane* and the *content pane*.

Note: You can tell the stage of a review by looking for the protocol or review icon at the top left of the *outline pane*.

When you undertake a Cochrane review, the protocol must be completed and submitted for publication before moving on to the review stage. For this tutorial, to convert the file from the protocol stage to review stage, use the following directions:

1. Click the **Properties** button  in the top-right of the *content pane*.
2. In the **Properties** window, select the **General** tab.
3. Note the automatically generated **Version Number** – this number will change each time the review is checked into Archie. Note that you do not need a separate **Review Number** unless your Review Group provides one.
4. Enter the **Version Description** ‘First draft’. When choosing a version description, make sure it is brief and clearly indicates what’s different about this draft, e.g. ‘Edited by [your initials]’ or ‘Methods section complete’.
5. For the **Stage**, select **Full Review**.
6. Click **OK**. Note that all of the section headings of the review are now activated for inclusion (there are no light grey section headings).

Saving a review

When you have finished working on a review from Archie, we recommend that you always save the review by checking it back into Archie using the Check In  button on the toolbar. This ensures the most recent version of your review is safely stored and can always be accessed from Archie by your co-authors and editors.

We recommend that you **do not** keep saved reviews on your computer, to avoid confusion with different versions. Be careful not to use a version of your review from your computer if a more recent version is available in Archie. If for some reason you do need to save a review to your computer, such as this tutorial review which is not from Archie, go to the **File** menu and select **Save As**, or click the **Save**  button on the toolbar. You can open a saved review by going to the **File** menu and selecting **Open**, or clicking the **Open** button  on the toolbar.

RevMan also creates automatic backups of your review. See the **Help** menu for information on retrieving a backup file.

Adding authors to a review


For a real review, your Review Group will add your name as an author of the review when they create the file, along with any known co-authors. You may wish to add additional authors later.

There are two ways to add authors to a review:

- **Add an author based on their record in Archie** – this option allows you to add authors who are already listed in the Archie database. You will need to have internet access and an Archie user account for this option.
- **Add a local record (name only) that can be linked to Archie later** – this option allows you to add authors that are not yet listed in Archie.

1. In the *content pane*, click the  icon next to the **Review information** heading to show further subheadings.
2. Click the  icon to open the **Authors** section. Note that three authors have already been added to this review.
3. In the *outline pane*, click the key icon  next to **Review information**.
4. Select the **Authors** heading.

5. Click the **Add Author**  button on the *outline pane* toolbar to open the **New Author Wizard**.
6. Select **Add a local record (name only) that can be linked to Archie later.**. Click **Next**.

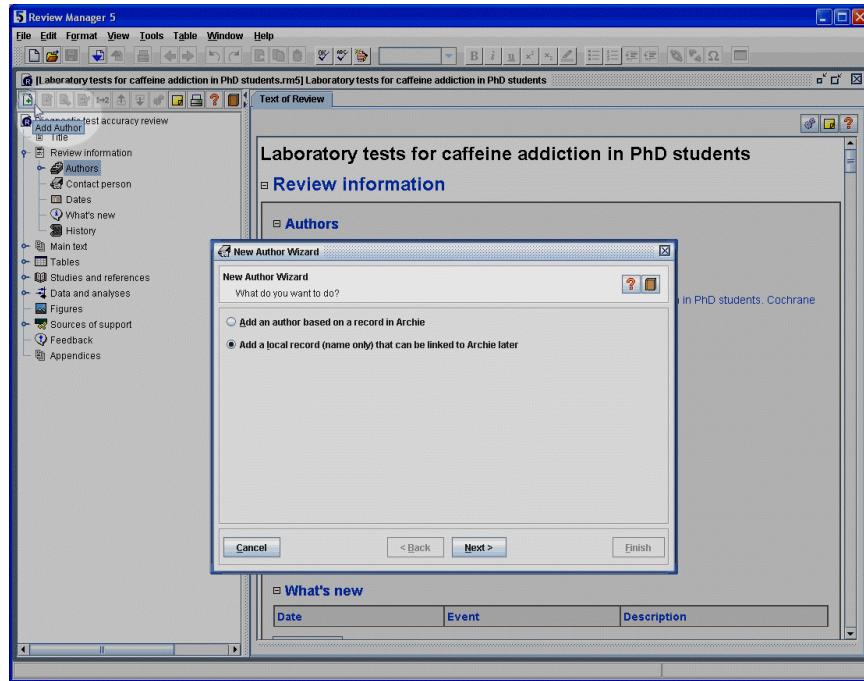

7. Enter your name in the **First Name** and **Family Name** fields.
8. Click **Finish**. You will now see your name listed next to the other authors of the review. Note that no affiliation is available for your name. Affiliations can only be added from records in Archie.

Under the author details in both the *outline* and *content panes* you will note that a **Contact person** can be specified for the review. By default, this is the first author, but another person (who may or may not be an author of the review) can be nominated by selecting the **Contact person** heading in the *outline pane* and clicking the **Edit Contact Person**  button on the *outline pane* toolbar. The contact person must be listed in Archie.

Note: To add new people to Archie or to make changes to contact details, please contact your Review Group Co-ordinator.

PART 2 – Working with text

Adding text to a review

Text can be copied and pasted from a word processor document, imported as an html file or typed directly into RevMan.

1. In the *content pane*, if it is not already expanded, click the  icon next to the **Background** heading.
2. Click the  icon next to the **Target condition being diagnosed** heading. A blank line will open under the heading.
3. Click to place the cursor in the blank line and type “Caffeine addiction is a common problem that affects productivity of PhD students, mainly due to too much talking at the coffee machine.”

Formatting text

RevMan text formatting options are available from the toolbar, including bold , italics , underline , superscript , subscript , and highlighting text  similar to most other word processing software. You cannot change the font or type size, as these are standard for all Cochrane reviews.


The **Insert Symbol** button  allows you to insert symbols into the text of the review, including in tables and author names. You can also insert a non-breaking space.

RevMan also allows bulleted  and numbered lists , including multi-level lists.

Making changes to headings

RevMan includes the standard headings for a complete Cochrane review, but new headings can be added and some standard headings can be deactivated.


1. In the *outline pane*, click the key icon  next to **Main text** to reveal the subheadings.
2. Click the key icon  next to **Background**.
3. Any heading with the  icon is optional and can be deactivated. Right-click the **Alternative test(s)** heading and select **Deactivate Heading** (alternatively, select the heading **Alternative test(s)** then click the  button on the *outline pane* toolbar and select **Deactivate Heading** from the drop down menu).
4. Note that the heading appears deactivated (light grey) in the *outline pane* and the text of the heading has been removed from the *content pane*.
5. In the *content pane*, click to place the cursor at the end of the text in the **Target condition being diagnosed** section and press **Enter** to create a new line.
6. On the toolbar, you will see a drop-down box showing the text style currently in use. By default, text is in the **Normal** style. Click the drop-down arrow and select **Heading 3**.
7. Type the text ‘Reference standard’. Note that formatting of Heading 3 is bold and italics.
8. Press **Enter** again. Note that the new blank line is in the **Normal** style again.

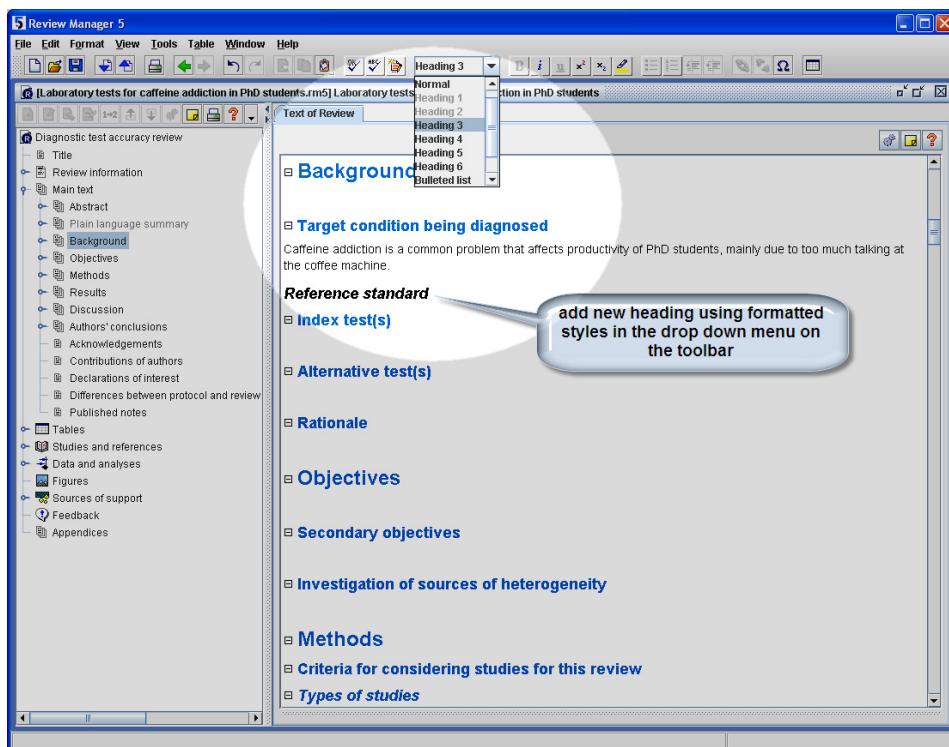

Adding notes

You can add notes to assist in drafting the review. These notes will not appear in the published version of the review.


1. In the *content pane*, make sure the cursor is still in the **Background** section.
2. Click the **Notes** button on the *outline pane* toolbar (NOT the **Notes** button in the *content pane*). This will open the **Notes** window.
3. Enter a note, 'This text needs revision.' You will see the user details you entered in the **Preferences** window at the beginning of the tutorial appears at the bottom of the **Notes** window.
4. Click the **Close** button to save the note and close the window.
5. A yellow **Note** icon appears next to the relevant heading in both the *content pane* and *outline pane*.
6. To see all the notes for a review while you are working, go to the **View** menu and select **Notes**. A separate window pane will appear below the *content pane* showing the notes for each section.
7. In the *content pane*, click to place the cursor back in the **Background** section. The notes area will automatically display the note you have entered in that section.
8. Click the icon at the top left of the notes area to minimise it, or go to the **View** menu and select **Notes** to close the notes pane.
9. In the *content pane*, click the **Note** icon to open the **Notes** window again.
10. Click **Delete** to delete the note.

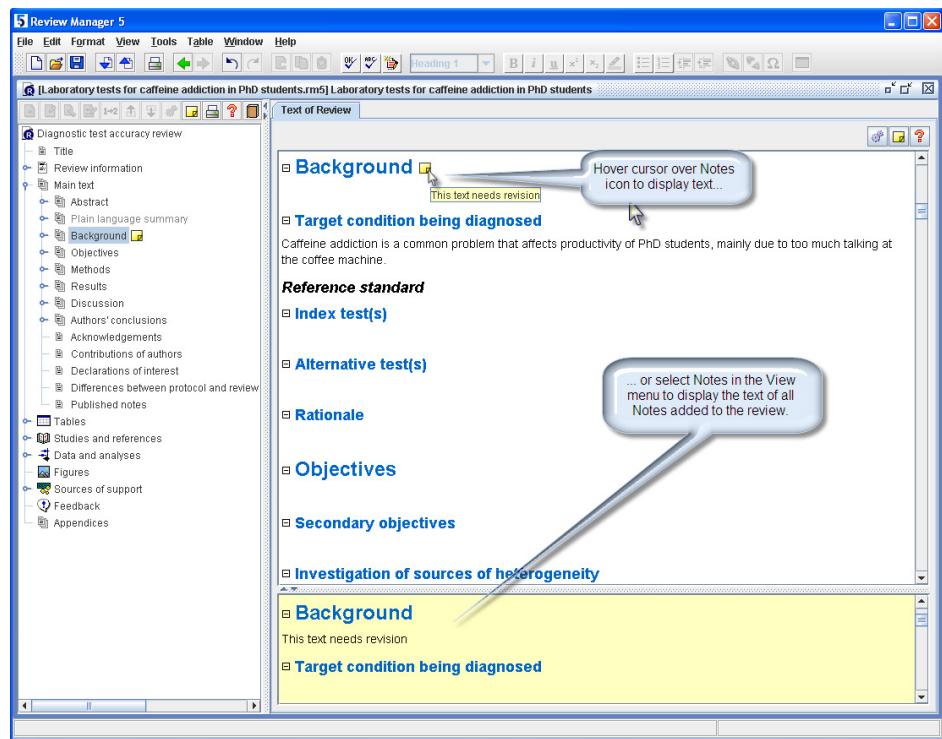

Tracking changes in the text

RevMan has a **Track Changes** system to help you manage editing in the text.

1. From the **Tools** menu, select **Track Changes**, or click the **Track Changes** button on the toolbar. Note that the **Track Changes** button is now activated.
2. In the *content pane*, go to the **Background** section and replace the word ‘problem’ with the word ‘condition’. The old text will appear in red with a line through it, and the new text will appear in green.
3. In the same section of the review, change the word ‘at’ to ‘by’. Leave the cursor in the text after your changes.
4. From the **Tools** menu, select **Accept or Reject Changes**. This will open the **Accept or Reject Changes** window.
5. In the window you just opened, click the **Find** button with the arrow pointing left to search backwards toward the beginning of your document.

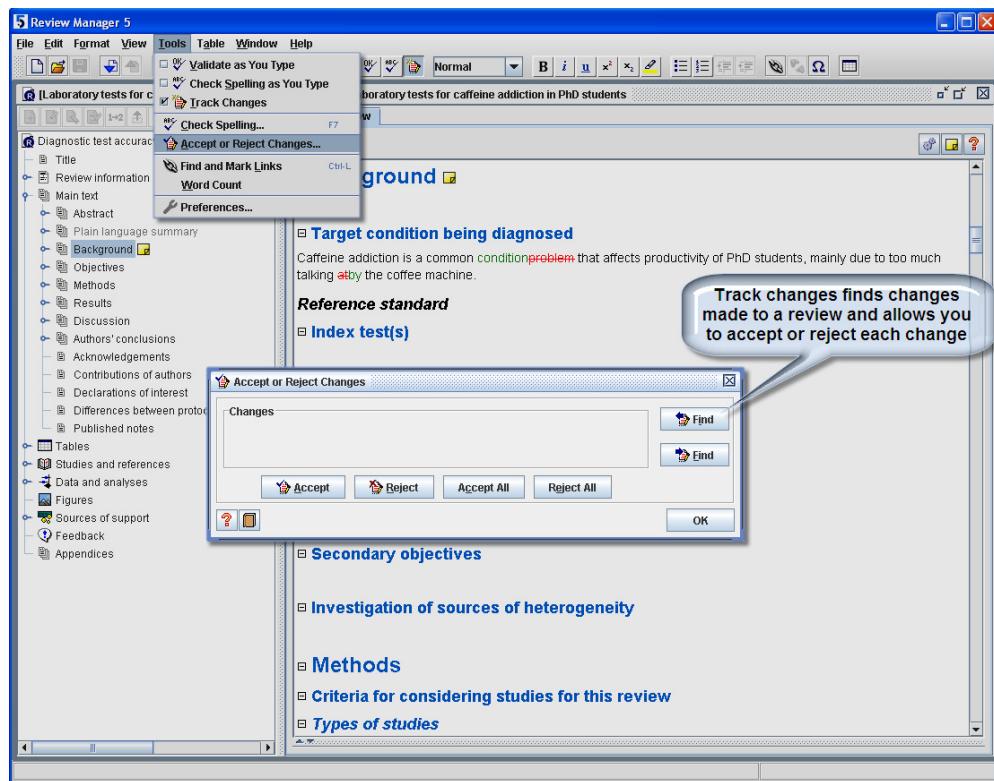

6. The first change found will be the word ‘function’. The word will be highlighted in the *content pane*. You may need to move the **Accept or Reject Changes** window out of the way to see the highlighted section of text. Click **Accept** to accept the change.
7. RevMan will automatically find the next change and highlight the relevant text. Click **Accept** again to accept the deletion of the word ‘at’.
8. Reject the rest of the changes you made by clicking the **Reject All** button.
9. Click **OK**.
10. In the *content pane*, change the word ‘problem’ to ‘condition’ again.
11. Hover the mouse pointer over the change. A popup message will appear showing the user details and when the change was made.
12. Right-click each change (both the added word and deleted word) and select **Accept Change**.
13. Turn off the **Track Changes** feature by clicking the **Track Changes** button on the toolbar.

Checking your spelling

RevMan can either check your spelling as you type, or check all the spelling in your review at the one time.

1. From the **Tools** menu, select **Preferences**.
2. Select the **Spell Checking** tab.
3. Select which dictionaries you would like RevMan to use, including **UK**, **US** or **Canadian English**. You can select more than one dictionary if you wish.
4. Click **OK**.
5. In the *content pane*, create a spelling mistake by deleting an ‘m’ from the word ‘common’.

6. Open the **Tools** menu and select the **Check Spelling** option. This will open the **Check Spelling** window.
7. Depending where your cursor was placed, RevMan may ask whether you wish to continue checking from the beginning of the document. Click **Yes**.
8. RevMan will highlight the spelling mistake and suggest ‘common’ as the correct spelling. Click **Change**.

Note: The **Check Spelling** tool cannot correct spelling errors in some sections of the review text, such as author's names and contact details or the **Data and Analyses** section.

9. A notice will appear to say that the spell check is complete. Click **OK**.
10. Click **OK** again to close the **Check Spelling** window.
11. To check spelling continuously while you type, go to the **Tools** menu and select **Check Spelling as You Type**, or click the **Check Spelling as You Type**  button on the toolbar.
12. Delete the ‘m’ from ‘common’ again, and click to place the cursor elsewhere in the text. The spelling mistake will be underlined in red.
13. To correct the word, right-click it and select the correct spelling from the options available.

Note: The **Check Spelling** tool identifies the "s" in constructions like "Index Test(s)" and "Comparison Test(s)" as spelling mistakes.

PART 3 – Adding studies and references

There are two kinds of references that can be added to a review:

- **References to studies** – these are references to studies that have been either included or excluded from your review. A study can have more than one reference, for example if there has been more than one published paper, or if you have received additional information from the study author.
- **Other references** – these are any other source you would like to cite in your review, such as background articles, methodological references, software, etc.

Adding studies to a review

When you are ready to add included and excluded studies to your review, you will need to create a Study ID for each study, and then add the associated references. For Cochrane reviews, the Study ID is usually the name of the lead author of the main published paper and the year of publication. If you have more than one reference with the same name and year, you can add a letter to the year to differentiate them (e.g. Egberts 2001a, Egberts 2001b).

1. In the *outline pane*, click the key icon  next to Studies and references.
2. Click the key icon  next to the **Reference to studies** section. You will see that there are already some studies included.
3. Select **Included studies** and click the **Add Study**  button on the *outline pane* toolbar. This will open the **New Study Wizard**.
4. Enter the **Study ID ‘Dorlas 2006’** and click **Next**.
5. From the **Data source** drop-down list, choose **Published and unpublished data**. Click **Next**.
6. The year the study was conducted will appear by default as ‘2006’. Click **Next**.
7. If your study has an ID number, such as a document ID, you can enter it here by clicking the **Add Identifier** button, selecting an ID type from the drop-down list, and entering the ID.
8. Click **Next**.
9. Select **Add another study in the same section** and click **Finish**.
10. Repeat steps 4-9 and add the following studies to the tutorial review:

Egberts 2001


Levett 2007

Smit 2006

11. When you have added the details for the last study, instead of step 9, select **Nothing** and click **Finish**.
12. In the *outline pane*, click the key symbol  next to **Included Studies** to see the new studies listed.

Adding references to a study

1. In the *outline pane*, select the **Dorlas 2006** study.
2. Click the **Add Reference**  button on the *outline pane* toolbar. This will open the **New Reference Wizard**.

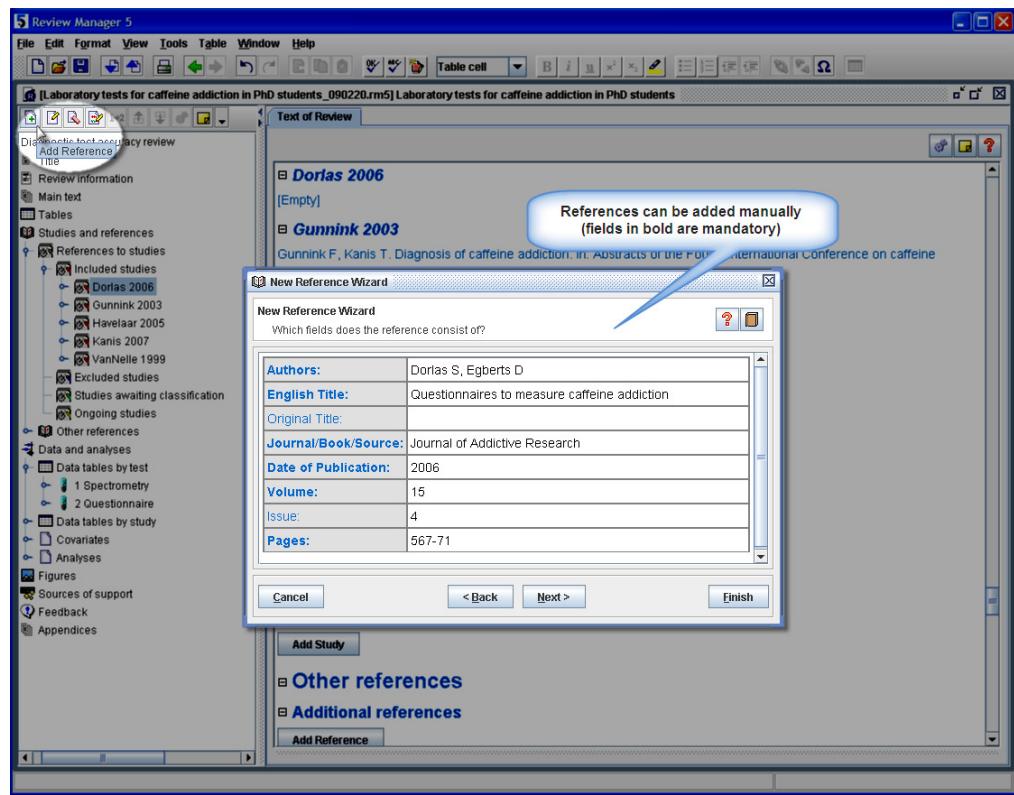

By default, the Reference Type is set to Journal article (click the drop-down list to see the other options available). Click Next.

3. Enter the following details in the form provided (note that fields appearing in bold are mandatory):

Authors: Dorlas S, Egberts D

English title: Questionnaires to measure caffeine addiction

Journal/book/source: Journal of Addictive Research

Date of publication: 2006

Volume: 15

Issue: 4

Pages: 567-71

Note: You can right-click in the **Journal/book/source** field and select **Choose from list** to select your journal from a list of titles.

4. Click **Next**. As for a study, you will have the option to add an ID number for the reference, such as a MEDLINE reference number.
5. Click **Finish**. The reference will now appear in the *content pane*. In the *outline pane*, click the key symbol next to **Dorlas 2006** to see the new reference listed.

Note: To add an **Other reference** to your review, select **Additional references** in the *outline pane* and follow steps 2 to 6 above. You will need to give your reference a Reference ID, which should have the same format as a Study ID (first author and year of publication).

Editing an existing reference

1. In the *outline pane*, double-click the reference listed under the **Dorlas 2006** study (or right-click and select **Edit Reference**). This will open a new tab in the *content pane* showing all the details you entered.
2. Open the **Reference Type** drop-down box and change the type from ‘Journal article’ to ‘Book’. Note that the fields have changed and some rows are now shaded. Shaded rows will not appear in your published review.
3. Change the reference information to the following (you need to scroll down to see some of these fields):

Authors: Dorlas S, Egberts D


Journal/book/source: Measuring caffeine addiction

Date of publication: 2006

Publisher Name: Grinder Press

City of publication: London

4. At the bottom of the reference tab, note that the reference is displayed as it will appear in the published review. Note that this display changes as you make the changes to the reference. Check that the reference appears correctly.

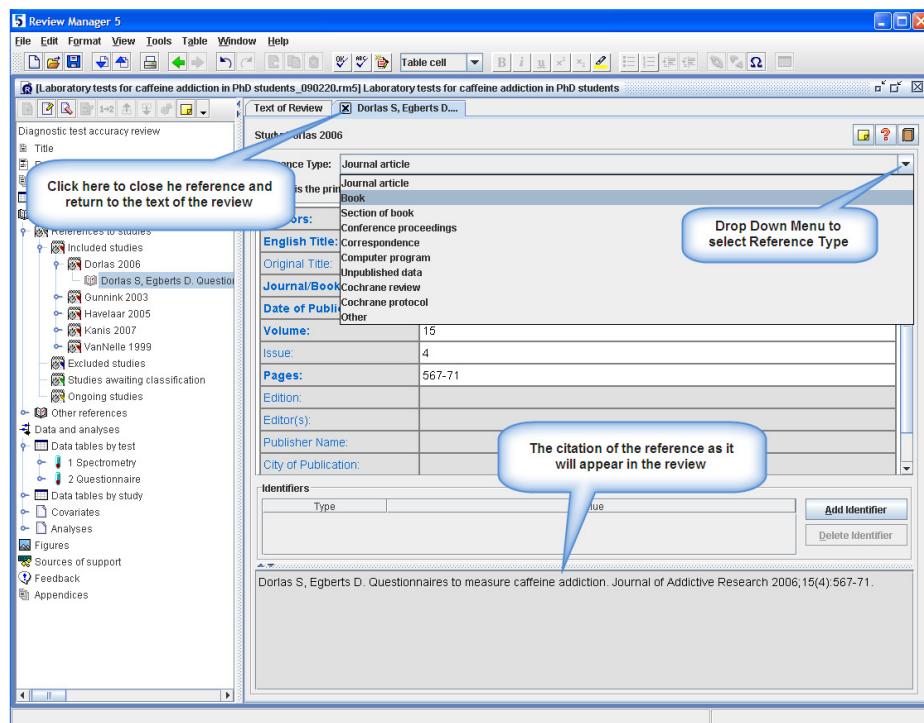

5. Go to the tab header and click the **X** button next to the first author’s name to close the reference tab (**do not** click the **X** button at the top right corner as this will close the review).


Note: For further details on how to enter different reference types, see the RevMan Help or the *Cochrane Style Guide*, available from the **Help** menu.

Importing references from a text file

If you have already entered your references in another program, you can avoid typing them in again by importing them from a text file. A sample text file has been provided with RevMan 5 for this tutorial.

Note: For details on how to generate a text file, see the help files in your reference management software. RevMan can read text files created from most major reference management software using the RIS format. You may need to check that you are using the right filter when exporting your references to a text file.

1. Go to the **File** menu and select **Import**, then **References**. This will open the **Import References Wizard**.
2. Click **Next**.
3. RevMan should automatically open the same folder we used to open the sample review: **c:\program files\review manager 5\doc\tutorial**. Select the file '**DTA Study References.txt**' and click **Open**.
4. You will see a preview of the text file. Click **Next**.
5. By default, **RevMan Format** is selected as the format of the text file (other options include RIS or PubMed MEDLINE formats). Click **Next**.
6. By default, RevMan will save references to the **Classification pending references** section, under **Other References** in the *outline pane*. References in this section will not appear in the published review (for this reason they are shown in grey). Click **Next**.
7. If you wish, you can assign different locations for each individual reference, or exclude specific references from the text file while importing the rest. For this tutorial, leave the settings as they are and click **Finish**.

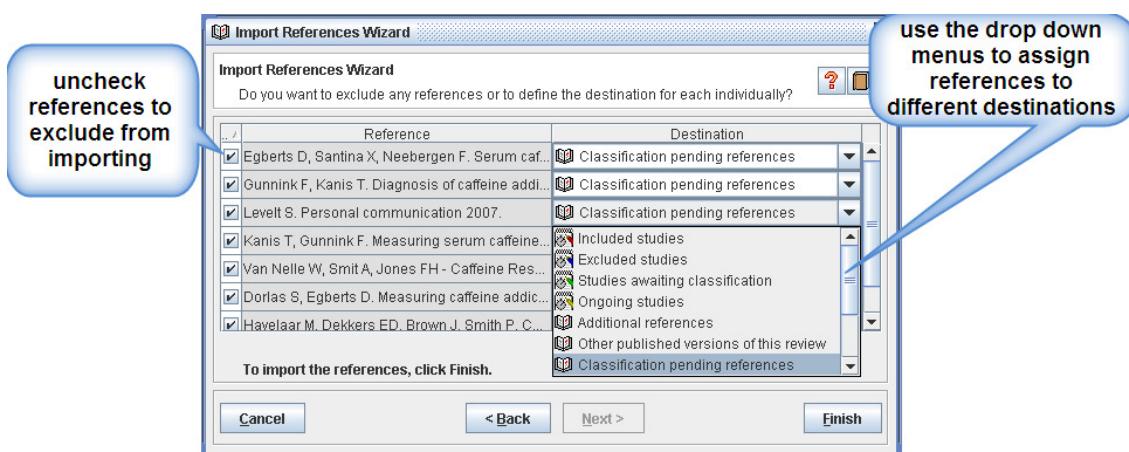

8. RevMan will report that 8 references have been imported. Click **OK**.
9. In the *outline pane*, use the key icons to open the **Other references** section and then the **Classification pending references** section. You will see the references you imported.

Moving references

There are several ways to move references to the appropriate study.


1. In the *outline pane*, in the **Classification pending references** section, select the reference **Egberts 2001**.

2. From the **Edit** menu, select **Cut** or click the **Cut** button  on the toolbar.
3. Use the key icons  to open the **Included studies** section, and open the **Egberts 2001** study. An [Empty] reference has been added here. To add a reference manually, you can edit the empty reference using the instructions in the previous section.
4. Select the **Egberts 2001** study.
5. From the **Edit** menu, select **Paste**, or click the **Paste** button  on the toolbar.
6. The reference will be deleted from the **Classification pending references** section and added to the **Egberts 2001** study. Note that the [Empty] reference has gone.
7. Go to **Classification pending references** and right-click **Gunnink 2003**.
8. Select **Move to**, then **Included studies**. As there was no study with a name matching the reference, RevMan will automatically create a new study in the **Included studies** section with the ID **Gunnink 2003**.
9. Go to the **Classification pending references** section and select **Kanis 2007**.
10. Click and drag the reference to the **Kanis 2007** study.
11. Repeat steps 1 to 6 and 10 to move the remaining references to the corresponding studies.
12. In the *outline pane*, click the key icons  to open the **Included studies** section, and then open the **Egberts 2001** study.
13. Select the first reference and double-click to open the reference tab.
14. In the *content pane*, check the **This is the primary reference for the study** box.
15. In the *outline pane*, note that the book icon beside the primary reference appears in black , and the icons beside the other references are grey .

Adding links to references in the text

Rather than typing citations in the text or including footnotes, RevMan allows you to include active hyperlinks to references in the text of your review.

1. In the *content pane*, scroll up to find the **Target condition being diagnosed** section.
2. Click to place the cursor after the word ‘machine’.
3. From the **Format** menu, select **Insert Link**, or click the **Insert Link**  button on the toolbar. This will open the **Insert Link** window.
4. From the list, select the study **Gunnink 2003**. Note that you can use this same method to link to included and excluded studies, figures, tables, sections of the review text and external web links.
5. Click **OK**. You will see a blue hyperlink to the reference appear in the text. This link will be active after the review is published in *The Cochrane Library*.
6. Add a space before the link and parentheses () before and after the link, so that it appears as “...the coffee machine ([Gunnink 2003](#)).”

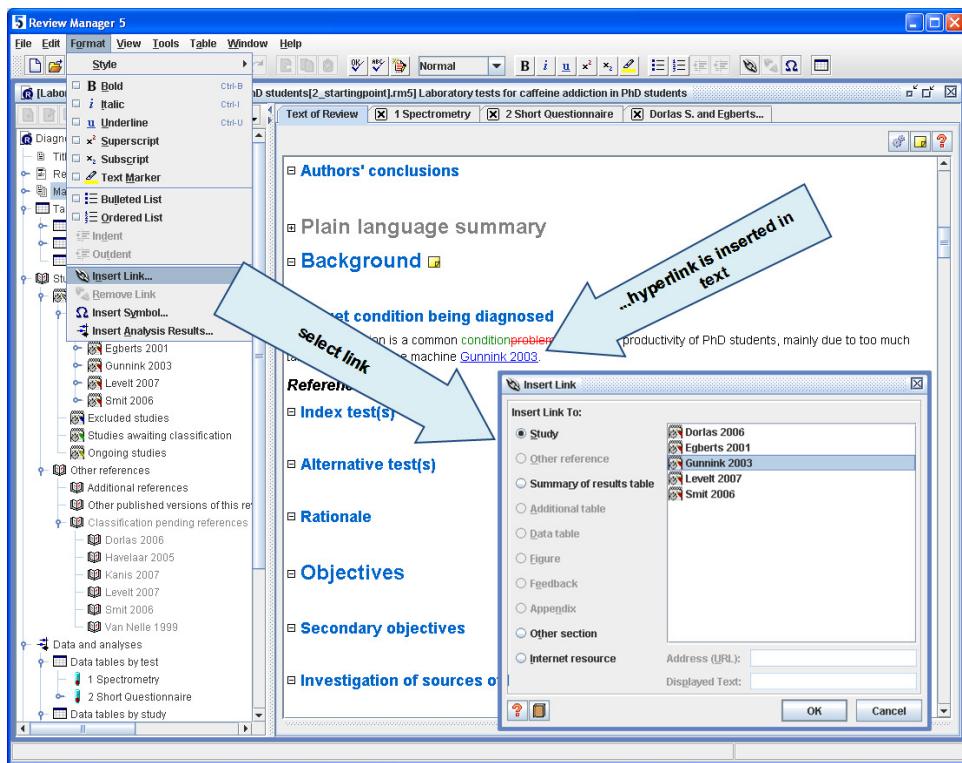

Removing a link

To remove a link you have already inserted in a review, you can either delete it as you would any other text, or use the **Remove Link** feature.

1. In the *content pane*, click to place the cursor in the link to **Gunnink 2003**.
2. Click the **Remove Link**  button on the toolbar, or right-click and select **Remove Link**. Note that the text is no longer blue.

Adding multiple links automatically

Instead of inserting one link at a time, RevMan can create all your links automatically by searching through the text to identify studies and other references.

1. In the *content pane*, select the sentence in the **Target condition being diagnosed** section, including the text '(Gunnink 2003)'.
2. From the **Tools** menu, select **Find and Mark Links**.
3. A notice will be displayed to let you know that one link has been created. Click **OK**. You will notice that the blue reference link has again been created.

Note: If you do not select a section of text before using the **Find and Mark Links** tool, RevMan will give you the option to check the whole review. Note that this tool will only search through text that is expanded (that is, visible in the *content pane*). Click the  icon next to a section heading to expand and see the text.

PART 4 - Tables

Now that you have included references to studies in your review, you can add a description of each study. This is done in the text of the review, and also in two standard tables: the **Characteristics of included studies** table and **Assessment of methodological quality** table.

Characteristics of included studies table

The details of each study, such as the methods, participants, target condition, and index and comparator tests should be described in the **Characteristics of included studies** table.

1. In the *outline pane*, use the key icon  to expand the **Tables** heading, then the **Characteristics of studies** heading and the **Characteristics of included studies** heading. You will see a list of the studies you have created.
2. Select the **Dorlas 2006** study.
3. In the *content pane*, you will see the same study shown. Click the  icon next to **Dorlas 2006** to show the **Characteristics of included studies** table.
4. Type or copy and paste the following text into the characteristics table:


Clinical features and settings	Primary care; sleep deprivation, loss of concentration, hyperactivity.
Participants	405 Dutch PhD students; no further information provided
Study design	Cross-sectional design; single group with equal suspicion of caffeine addiction; consecutive enrolments; every student underwent all tests.
Target condition and reference standard(s)	Caffeine addiction, measured by coffee consumption in liters per day
Index and comparator tests	Index test is spectrometry; comparator is questionnaire.
Follow-up	Two weeks follow up. Three students withdrew; reasons not explained.
Notes	Sources of funding not stated.

Assessment of methodological quality table

Assessment of methodological quality tables are used to describe the methods used in each study, and whether the results are at risk of bias. For information on assessing the methodological quality and its relation with risk of bias, see the *Cochrane Handbook for Systematic Reviews of Diagnostic Test Accuracy*, available from the **Help** menu or <http://srdta.cochrane.org/en/authors.html>.

1. In the *content pane*, under the **Characteristics of included studies** table, find the **Assessment of methodological quality table** heading (you may need to click the  icon).
2. Click the **Properties**  icon next to the **Assessment of methodological quality table** heading.
3. You will see a list of items that should be included in an **Assessment of methodological quality** table. All eleven listed items have a tick symbol  next to it, indicating it is included in the table by default.

4. Only in very specific cases (see the Handbook for this), authors may decide not to include one of the items in the table. Deactivate **Relevant clinical information** by selecting it from the list and by clicking on the **Activated** box, which then becomes deactivated.
5. Click ‘OK’.
6. The **Relevant clinical information** item will now be removed from the table for every study in your review. (Please note, that you will lose any text that had already been entered in the description column!)
7. Set each of the drop-down boxes in the table to **Yes**, **No** or **Unclear**. Text can now be typed directly into the table, to explain why you chose Yes, No or Unclear. For now, leave the **Description** column blank.

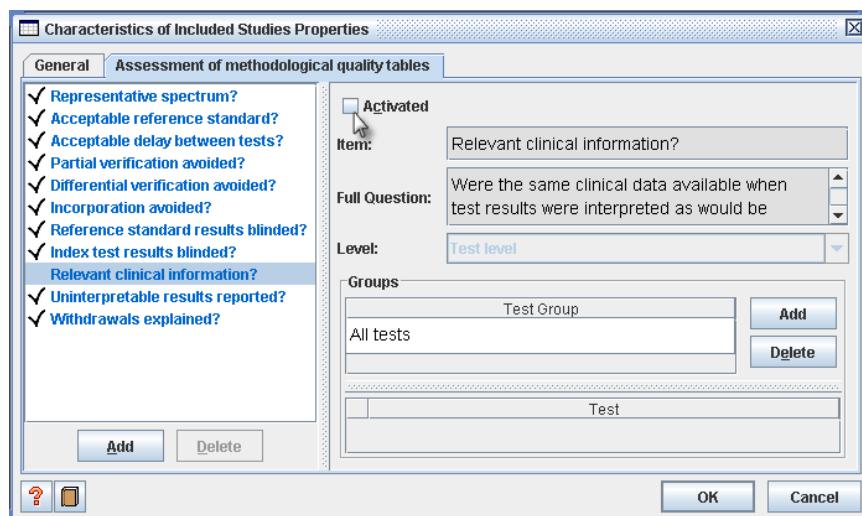


Some items are assessed for the study as a whole (e.g. **Acceptable reference standard**) while others must be assessed at the test level (e.g. **Acceptable delay between tests**). RevMan allows you to add separate information for each test or group of tests for each quality item. We will come back to that later, after we have added tests to the review (see Part 5, page 23).

Adding self-defined quality items

In some cases, it may be useful to add user-defined quality-items. For example, if you are assessing a commercially available device, it may be useful to check for sponsoring.

1. In the *content pane*, under the **Characteristics of included studies** table, find the **Assessment of methodological quality table** heading (you may need to click the  icon).
2. Click the **Properties**  icon next to the **Assessment of methodological quality table** heading.
3. Click on the **Add** button under the list of quality items. A new item will appear, called New Item, with a tick mark in front of it.

4. In the right-hand half of the window, in the **Item** field, change the name of the new item to '**Sponsoring precluded?**'.
5. In the **Full Question** field, enter the full question: 'Were none of the authors paid by the manufacturer and were none of the tests provided by the manufacturer?'. Whenever you add a new item to the table, make sure you phrase the question so that answering "Yes" indicates higher quality.
6. Click '**OK**'.
7. The **Assessment of methodological quality table** for each study will now include a row for the new heading.
8. Set the new drop-down box to **Yes, No or Unclear**. Text can now be typed directly into the table, to explain why you choose Yes, No or Unclear.

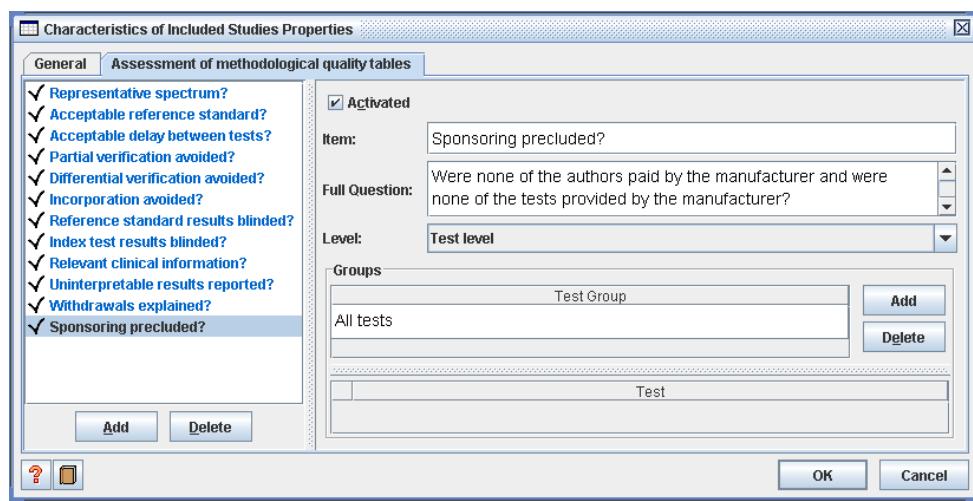

Note: In the **Assessment of methodological quality** tables, when you move your mouse over the displayed items, you will see the full question for this item appearing.

Item	Judgement	Description
Representative spectrum?	Unclear	
Acceptable reference standard?	Yes	Was the spectrum of patients representative of the patients who will receive the test in practice?
Acceptable delay between tests? (Spectrometry)	Yes	
Acceptable delay between tests? (Questionnaire)	Yes	
Partial verification avoided?	Unclear	
Differential verification avoided?	Yes	
Incorporation avoided?	Yes	
Reference standard results blinded?	Unclear	
Index test results blinded?	Unclear	
Relevant clinical information?	No	
Uninterpretable results reported?	Yes	
Withdrawals explained?	No	
Sponsoring precluded?	Unclear	

Characteristics of excluded studies

Studies identified in your search that met the exclusion criteria of the review may be listed under **Characteristics of excluded studies**, along with the reason each study

was excluded. You do not need to list every study identified by your search, but it is useful to list those which might appear to readers to meet your inclusion criteria.

1. In the *outline pane*, select **Characteristics of excluded studies**.
2. Click the **Add Study**  button on the *outline pane* toolbar and enter the new study 'Bean 1988'. Click **Next**.
3. Select **Published data only**, and click **Next**.
4. Accept '1988' as the year the study was conducted, and click **Next**.
5. Click **Next** without adding any identifiers.
6. Accept '**Nothing**' as your next action, and click **Finish**. In the *content pane*, you will see the Excluded studies section and the study you have created.
7. In the *outline pane*, click the key icon  next to **Characteristics of excluded studies**.
8. Select the new study **Bean 1988**. In the *content pane*, you will see a table under this study with the heading **Reason for exclusion**.
9. In the table, enter the text, 'Coffee machine only contained decaf coffee.'

Additional tables


Additional information can be included in table format either in the text of the review, in an **Appendix** or as an **Additional table**.


To add a table into the text of the review:

1. In the *content pane*, scroll up to the **Background** section and place the cursor in the **Target condition being diagnosed** section.
2. Click the **Insert Table** button  on the toolbar. A 3x3 table will be created.
3. Right-click the top left cell. Note the table edit options.
4. Select **Toggle Heading/Cell**. The selected cell will be shaded, and any text entered will appear in bold. Text can now be entered directly into the table.

Adding Footnotes to tables

In the content pane, below each tables section, there is a place to add footnotes for these tables. To add Footnotes, click on the white space below the Caption heading and type in some text.

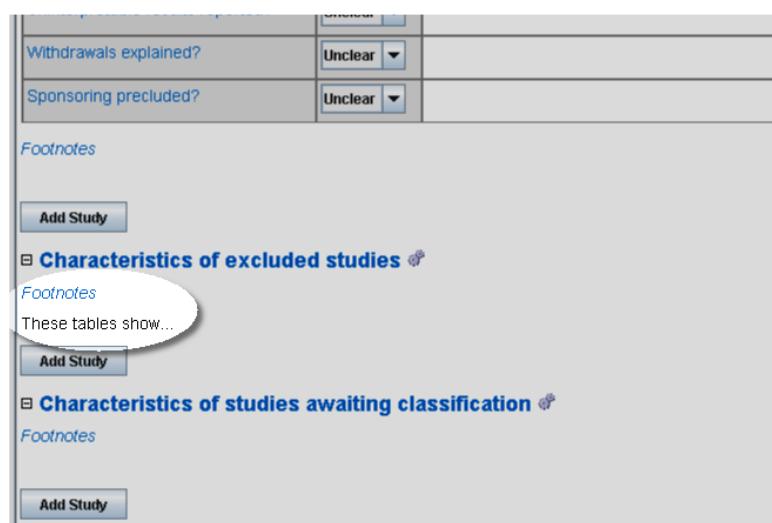

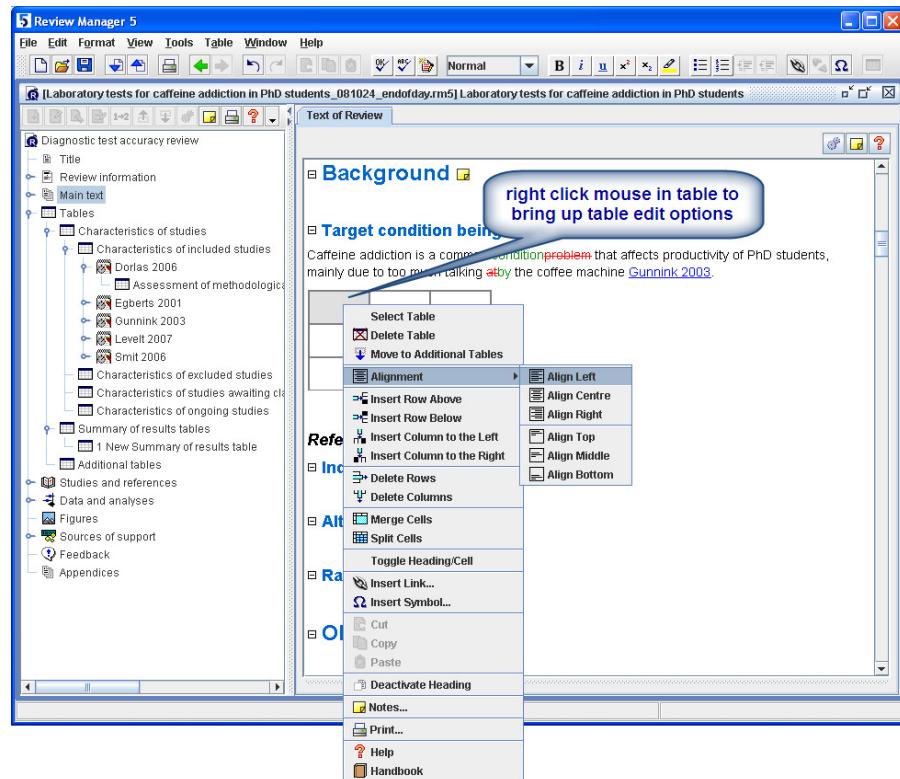

To add an **Additional table** that will appear at the end of the published review:

1. In the *outline pane*, under the **Tables** heading, select **Additional tables**.
2. Click the **Add Table** button  on the *outline pane* toolbar. This will open the **New Additional Table Wizard**.
3. Enter the title, ‘Decaffeinated coffee: chemical information’.
4. Click **Next**.
5. Accept the default number of rows and columns and click **Finish**. In the *content pane*, you will see the table in the **Additional tables** section. You can now enter text directly into the table.
6. In the *outline pane*, click the key icon  to see the new table listed under **Additional tables**.

PART 5 - Data and analyses


Once you have described your included studies, you can begin analysing the results. Before analysing the results, you will need to enter data and, if applicable, covariates. Entering data can be done in two ways. You can enter them by test, under **Data tables by test**. Using this option, you enter the data for one test at a time as reported in each study. Alternatively, you can enter the data by study, under **Data tables by study**. Using this option, you enter the data for one study at a time, for all tests reported in that study. This may be advantageous when the included studies analysed many tests.

When you have added included studies in your RevMan file, you will see the key icon  next to **Data tables by study**. That means that RevMan has listed all your included studies under that heading.

Adding tests

The first step in analysing diagnostic accuracy data in RevMan is to enter the tests that you are interested in. Reviews may include one or more tests.


1. In *outline pane*, select **Data tables by test**.
2. Click the **Add Test**  button on the *outline pane* toolbar. This will open the **New Test Wizard**.
3. Enter the **Name** ‘Spectrometry’ and, if applicable, the Full Name of the test (for example, if you have different manufacturers for spectrometers). Then click **Next**.
4. You can add a description of the test in the next pane. Click **Finish**. This will open an empty Data Tables tab.
5. Repeat steps 1-4 to add another test called ‘Questionnaire’.

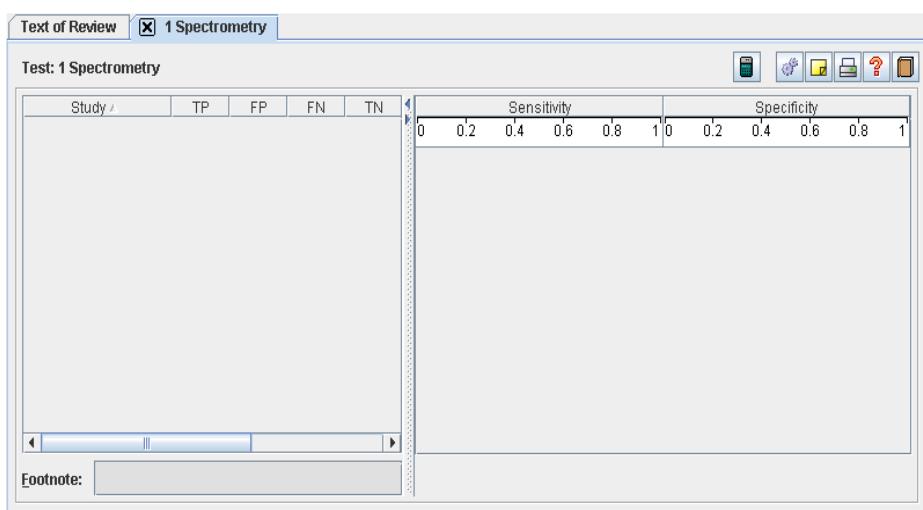

Entering data by test

When you are interested in only one or a few tests, or when you expect the included studies to report the results of only one or a few tests, it will probably be easiest to enter your data by test.


1. In the *outline pane*, click the key icon  to see the new test(s) listed under **Data tables by test**.
2. Select **Spectrometry** and click the **Add Test Data**  button on the *outline pane* toolbar. This will open the **New Test Data Wizard**.
3. You will now select the studies that report data for Spectrometry. By holding down the Ctrl/Command-button on your computer, you can select multiple studies at once.
4. Select the following studies:
 - Gunnink 2003
 - Havelaar 2005
 - Kanis 2007
 - Smit 2006
 - VanNelle 1999

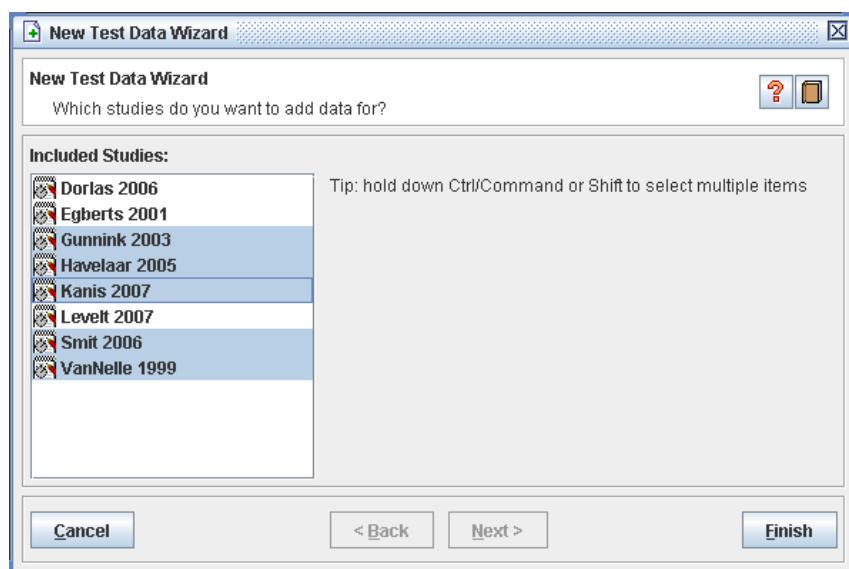

5. Click **Finish**, this will open the Data Tables ready for you to enter the selected studies listed. Beside each study is a row of cells ready for you to enter the numbers of true positive results (TP), false positive results (FP), false negative results (FN) and true negative results (TN).
6. Enter the following data into the table:

Study	TP	FP	FN	TN
Gunnink 2003	18	4	32	46
Havelaar 2005	22	17	1	38
Kanis 2007	60	87	16	371
Smit 2006	8	8	4	59
VanNelle 1999	9	2	5	65

7. RevMan will automatically calculate the sensitivity and specificity (with 95% confidence intervals) for each study. The results are also displayed in a forest plot. The sensitivity and specificity for each study are shown as a blue square, and with a horizontal line showing the confidence intervals.

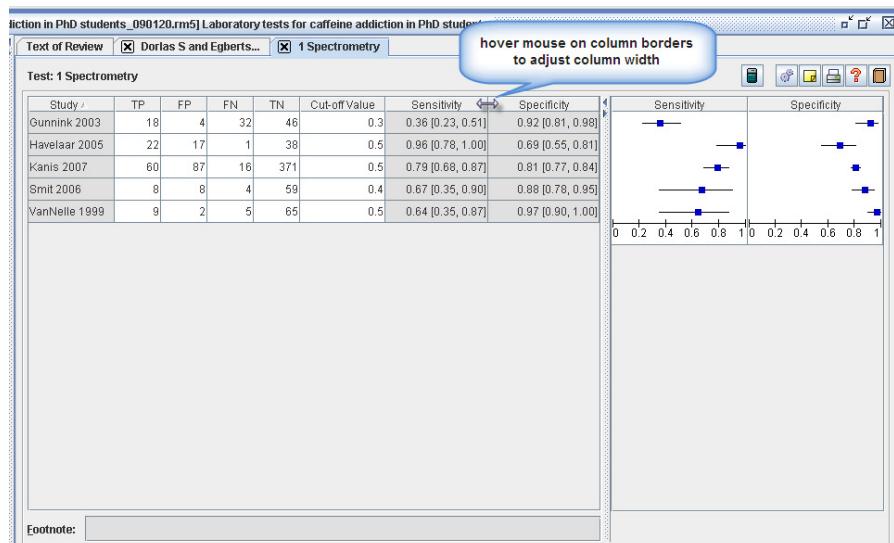

Note: For those of you who are familiar with intervention reviews, note the following. There is no pooled result shown for sensitivity and specificity, RevMan does not calculate weights for the forest plots, and no measures of heterogeneity are given. For more details on why these features are not provided, please refer to the *Cochrane Handbook for Systematic Reviews of Diagnostic Test Accuracy*, available from the Help Menu or at <http://srtda.cochrane.org>.


8. At the bottom of the screen under the results table is space to enter footnotes. In the results table, select the **Kanis 2007** study. In the **Footnote** field, enter the text 'Unpublished data provided by author'.
9. At the top of the *content pane*, click the **Text of review** tab. Note that RevMan has added a table with a brief summary of the results for this outcome, including the number of studies and participants.

Entering data using the calculator

Another way to enter your data is to use the calculator. This is particularly useful when you do not have the numbers of true and false positives and negatives, but you do have sensitivity, specificity and the total numbers of people with and without the target disease.

1. Select under **Data tables by test** the test called **Questionnaire** and click the **Add Test Data** button on the *outline pane* toolbar. This will open the **New Test Data Wizard**.
2. Select the following studies that report data for Questionnaire and click **Finish**:
 - Dorlas 2006
 - Egberts 2001
 - Gunnink 2003
 - Havelaar 2005
 - Levelt 2007
3. Select the first cell in the first row (Dorlas 2006) and click on the calculator icon . A new window opens.
4. Enter 0.75 in the Sensitivity box. Note that RevMan will only recognize numbers between 0 and 1, so you cannot enter 75%.

5. Enter 8 in the D+ box. D+ stands for ‘Diseased’. This means that there were in total eight people with caffeine addiction included in the study.
6. RevMan will now automatically calculate TP and FN.
7. Enter 0.67 in the Specificity box and 91 in the D- box. D- stands for ‘Not Diseased’. This means that there were in total eight people without caffeine addiction included in the study.
8. RevMan will now automatically calculate FP, TN and all other measures, including positive and negative likelihood ratios. It also checks for correctness. If the numbers seem to be correct (i.e. if the given parameters can be used to calculate round numbers of TP, FP, FN, TN), the boxes appear green. Otherwise the sensitivity box and/or specificity box will appear red.
9. Click ‘OK’. The cells in the first row will be filled with numbers and a forest plot of sensitivity and specificity will be displayed.

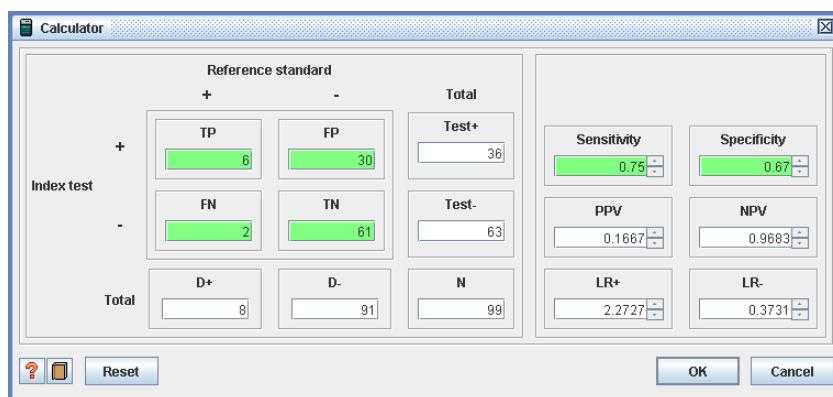


Entering data by study

When you expect the included studies to report the results of many tests at one time, you may want to enter your data tables by study.

1. In the *outline pane*, click the key icon to see the studies listed under **Data tables by study**.
2. Click the key icon next to **Egberts 2001** to see the tests that are already added to this study. Double-click on **Egberts 2001** to open the data tables for this study.
3. Enter these data into the table:

Study	TP	FP	FN	TN
Egberts 2001	19	9	6	51

4. In the *outline pane*, click the key icon next to **Gunnink 2003** to see the tests that are already added to this study. Double-click on **Gunnink 2003** to open the data tables for this study.
5. You will see that RevMan has already entered the data for Spectrometry, taking the data you entered under **Data tables by test**.

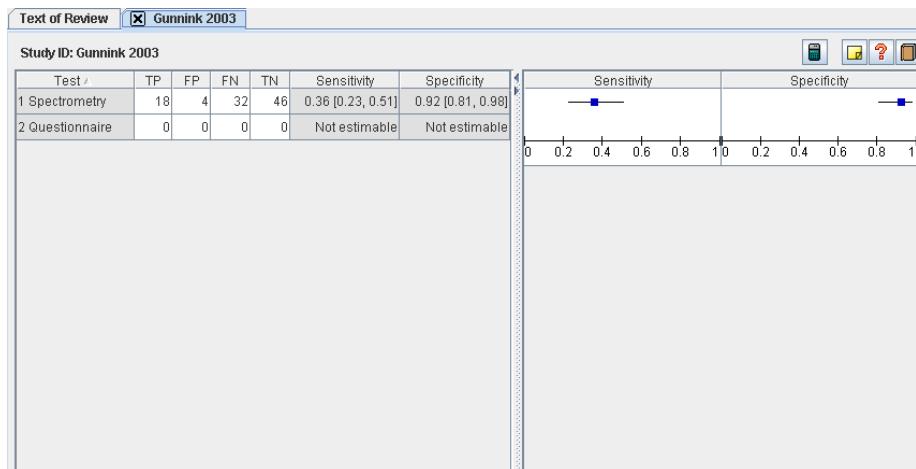

6. Enter the data for the Questionnaire:

Study	TP	FP	FN	TN
Gunnink 2003	10	15	2	53
Havelaar 2005	74	54	26	246
Levett 2007	9	23	1	51

Note: You can also enter the data by copying these from a spreadsheet. In that case, you select all the cells containing TP, FP, FN, TN and copy those cells. Then you right-click on the first left-hand cell in the RevMan data table and choose **Paste**. Make sure that the rows and columns in your spreadsheet and in RevMan follow the same order!

Note: When you update a value in the data tables, make sure that you press Enter/Return afterwards. If you do not press Enter, RevMan will not regard the number as really being entered and will not update the matching forest plots.

Assessment of methodological quality for two tests

Some quality items are assessed for the study as a whole (e.g. **Acceptable reference standard**) while others must be assessed at the test level (e.g. **Acceptable delay between tests**). RevMan allows you to add separate information for each test or group of tests for each quality item.

1. In the *content pane*, select the **Assessment of methodological quality table** heading.
2. Click the **Properties** icon next to the **Assessment of methodological quality table** heading.
3. By default ‘All tests’ is displayed in the **Test Group** area.
4. To assess the item **Acceptable delay between tests** separately for both index tests included, select this item and click the **Add** button next to the **Test Group** area.
5. Delete the default text in the upper box in the **Test Group** list, ‘All tests’, and enter ‘Spectrometry’. Change the default text in the lower box into ‘Questionnaire’.

6. The **Activated** box should still be checked.
7. Select **Spectrometry** in the **Test Group** list and check the **Spectrometry** box in the **Test** list.
8. Select **Questionnaire** in the **Test Group** list and check the **Questionnaire** box in the **Test** list.
9. See if the **Activated** box is still checked.
10. Click **OK**.
11. The **Assessment of methodological quality table** will now include a row for each activated heading, including separate rows on **Acceptable delay** for each test.
12. Set each of the drop-down boxes in the table to **Yes**, **No** or **Unclear**. Text can now be typed directly into the table, to explain why you choose Yes, No or Unclear.

□ [Assessment of methodological quality table](#)

Item	Judgement	Description
Representative spectrum?	Unclear ▾	
Acceptable reference standard?	Yes ▾	
Acceptable delay between tests? (Spectrometry)	Yes ▾	
Acceptable delay between tests? (Questionnaire)	Yes ▾	
Partial verification avoided?	Yes Yes Unclear No	
Differential verification avoided?	Yes ▾	
Incorporation avoided?	Yes ▾	
Reference standard results blinded?	Unclear ▾	
Index test results blinded?	Unclear ▾	
Relevant clinical information?	No ▾	
Uninterpretable results reported?	Yes ▾	
Withdrawals explained?	No ▾	

Covariates

It is possible to add covariates as well, for example, the cut off values of the tests that are being assessed.

1. In *outline pane*, right-click **Covariates** and choose **Add Covariate**. This will open the **New Covariate Wizard**.
2. Enter the **Name** ‘Cut-off Value’ and click **Next**.
3. Select **Test Level** and click **Next**.
4. There are two types of covariates: **Categorical** or **Continuous**.
5. Select **Continuous**. For continuous covariates, you can leave the **Default Value** blank or you can enter, for example, 0 for default.
6. Click **Finish**.
7. In the **Data Tables** (for both tests), a column will appear for **Cut-off Value**.
8. If you left the **Default Value** blank, you will see an empty column, otherwise, this column will be filled with your default value.
9. Enter the following values in the tables (either via **Data tables by test** or **Data tables by study**).


Study	Test	Cut-off Value
Gunnink 2003	Spectrometry	0.3
Havelaar 2005	Spectrometry	0.5
Kanis 2007	Spectrometry	0.5
Smit 2006	Spectrometry	0.4
VanNelle 1999	Spectrometry	0.5

10. Note that both Covariates appear in the Data Table of all included Tests. We leave cut-off value for questionnaires blank.
11. For the questionnaire, we will add another covariate, called ‘Questionnaire Type’.
12. In *outline pane*, right-click **Covariates** and select **Add Covariate** again.
13. Enter the **Name** ‘Questionnaire Type’ and click **Next**.
14. Select **Test Level** and click **Next**.
15. Select **Categorical** and click **Next**.
16. For categorical covariates, you must define the different categories. Click **Add** to create a category.
17. Delete ‘**New category**’ and replace it with ‘**Older List**’.
18. Click **Add** to create the second category. Replace the name ‘**New Category**’ with ‘**Newer List**’.



19. Click **Finish**.
20. In the **Data Tables** (for both tests), a column will appear for **Questionnaire Type**. In the **Questionnaire** data table, use the drop down menus to assign the questionnaire types as follows:

Study	Test	Questionnaire Type
Dorlas 2006	Questionnaire	Older List
Egberts 2001	Questionnaire	Newer List
Gunnink 2003	Questionnaire	Newer List
Havelaar 2005	Questionnaire	Older List
Leveld 2007	Questionnaire	Newer List

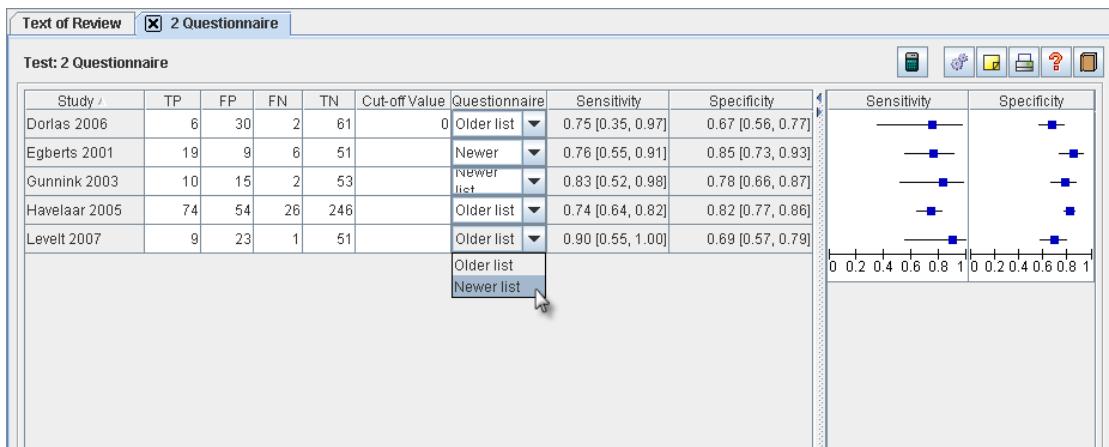


Note: Although you may wish to explore the effect of certain quality items, these need not to be entered as covariates. The quality items are already available as part of your Assessment of methodological quality tables.


Analyses

RevMan is not able to perform meta-analysis or calculate summary estimates for DTA studies. You will need external software to do this. RevMan is capable of graphical analysis of the data at hand, either on its own or with the input of parameters from external software.

1. In *outline pane*, select **Analyses**.
2. Click the **Add Analysis** button on the *outline pane* toolbar. This will open the **New Analysis Wizard**.
3. Enter the Name ‘Spectrometry’ and click **Next**.
4. Select **Single test analysis**, and then select **Spectrometry** from the list of available tests. Click **Next**.


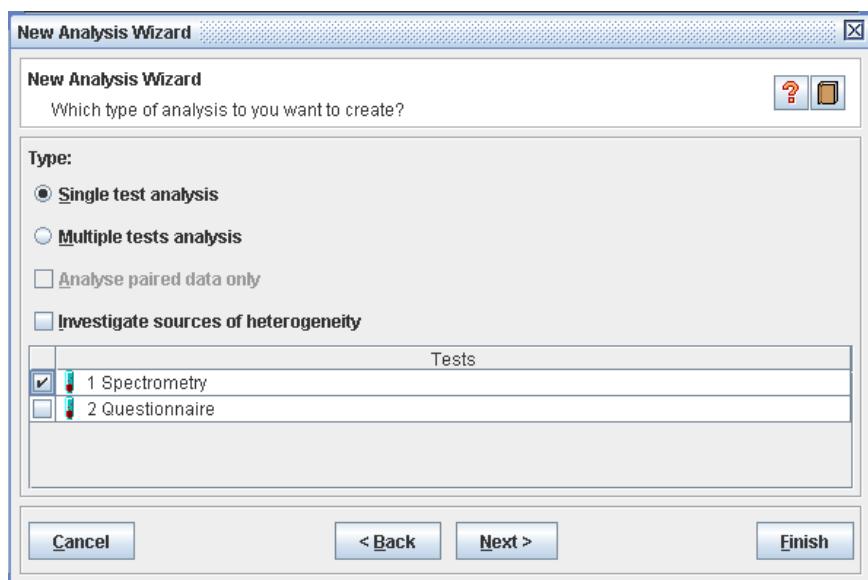


5. Click **Finish**.
6. An Analysis pane will open, showing three different areas:
 - a. The **Forest plot area**
 - b. The **SROC plot area**
 - c. The **Externally Calculated Parameters area**.

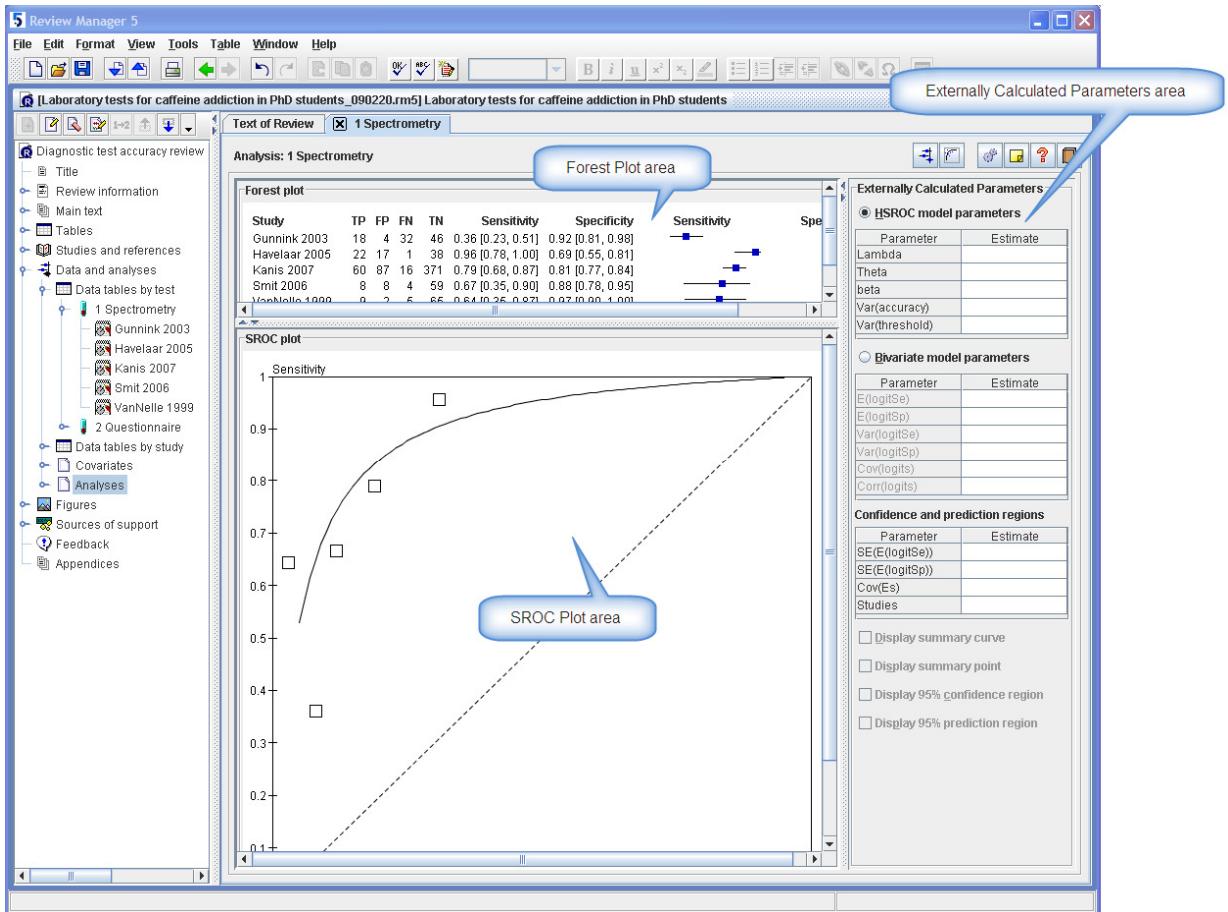


a. Changing Forest Plot properties

1. In the outline pane, click the key icon to see the analyses listed under **Analyses**.
2. Click the **Properties** button in the top right corner. This will open the **Analysis Properties** window.
3. Click the **Forest Plot** tab.
4. Select **Cut-off Value** as the **Covariate Displayed on Forest Plot**. Click **Apply** and then **OK**. Note that a column is added to the forest plot.

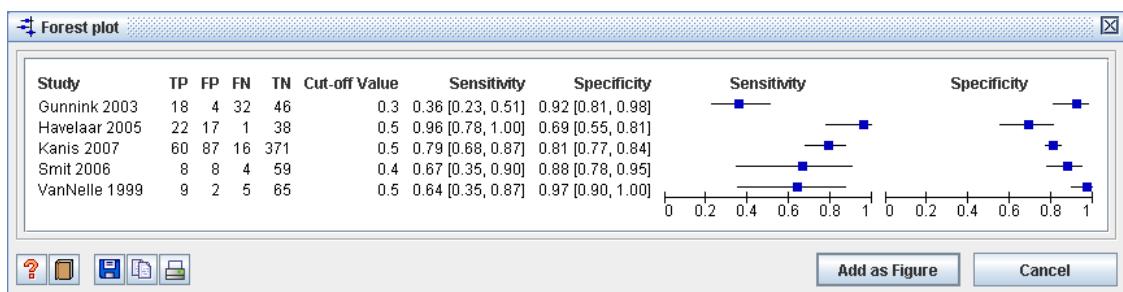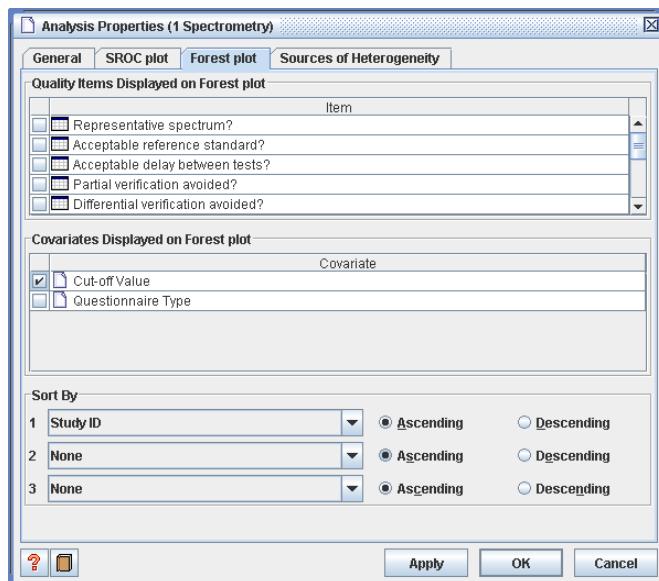


b. Changing SROC Graph properties

1. Click the **Properties** button in the top right corner to open the **Analysis Properties** window again.
2. Click the **SROC Plot** tab.

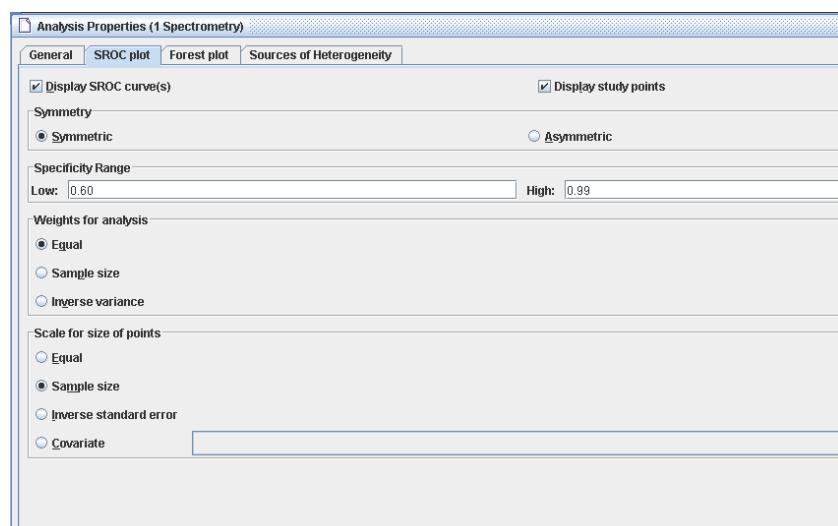

3. Select **Study Size** as the **Scale for size of points**. Click **Apply**. Note that the forms and size of the points in the graph have changed.
4. Under **Symmetry**, select **Asymmetric**. Click **Apply**. Note that the shape of the curve changed.
5. Set the low end of the **Specificity Range** at 0.60 and the high end at 0.99. Click **Apply**. Note that the SROC curve covers a different range in the plot.
6. Click **OK**.

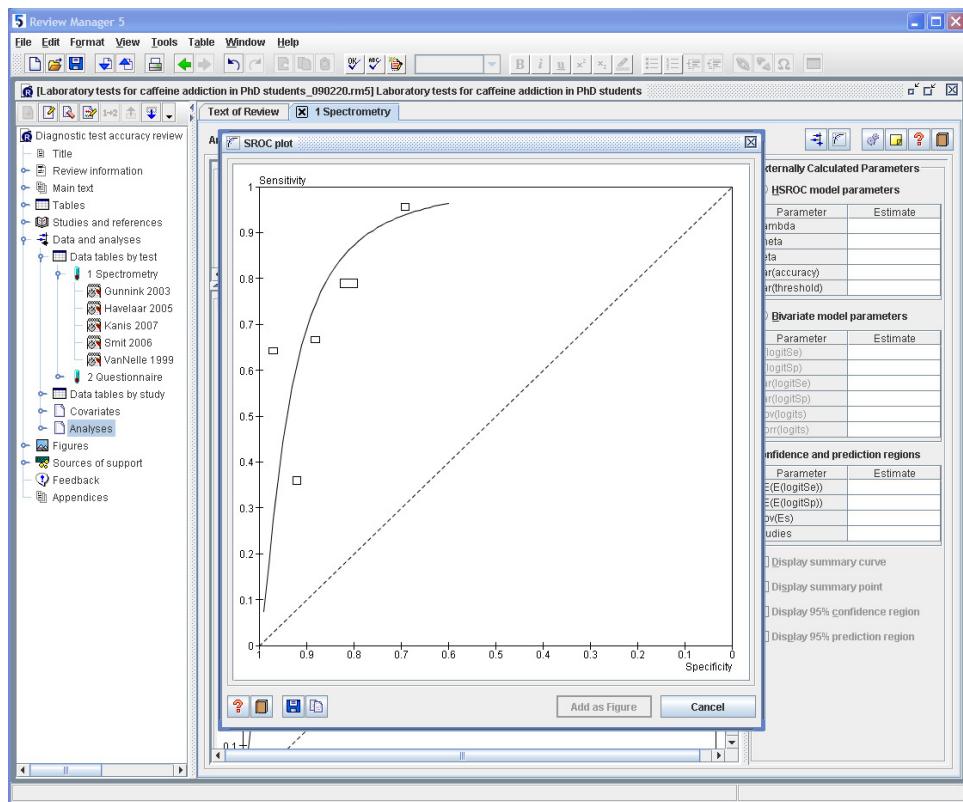

c. Entering Externally Calculated Parameters


To obtain valid pooled estimates of sensitivity and specificity, and when statistical testing is needed (e.g. to test whether the summary ROC curves of two tests significantly differ), external software has to be used. The resulting parameters can be entered here. Please consult the DTA Handbook or a statistician to explain how these parameters should be used.

Test Comparisons

Most review questions will be comparative rather than looking at one specific test.


1. In *outline pane*, select **Analyses**.
2. Click the **Add Analysis**  button on the *outline pane* toolbar. This will open the **New Analysis Wizard**.
3. Enter the **Name** ‘Comparison’ and click **Next**.
4. Select **Multiple tests analysis**.
5. Select both **Spectrometry** and **Questionnaire** and click **Next**.

6. Click **Finish**.

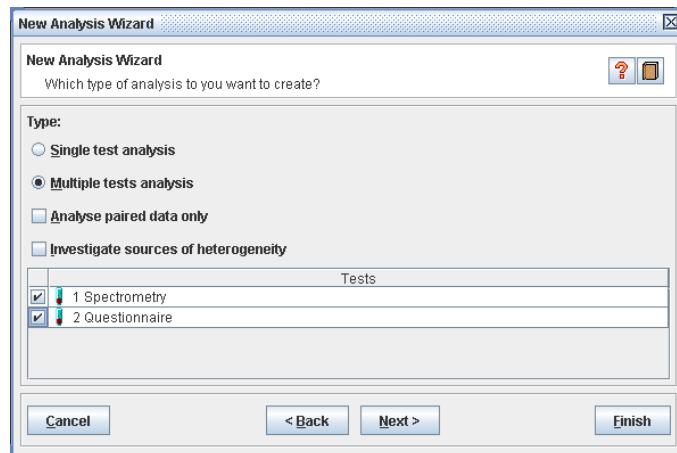

7. After clicking finish, the Analysis pane occurs and again you will see three different areas:


- The **Forest plot** area: scroll down to see two forest plots, one for each test.
- The **SROC plot** area: two different symbols with two different colours are now used. The black squares indicate the spectrometry studies and the red diamonds indicate the questionnaire studies.
- The **Externally Calculated Parameters** area: a drop-down menu appears to enable you to enter the externally calculated parameters for each test separately.

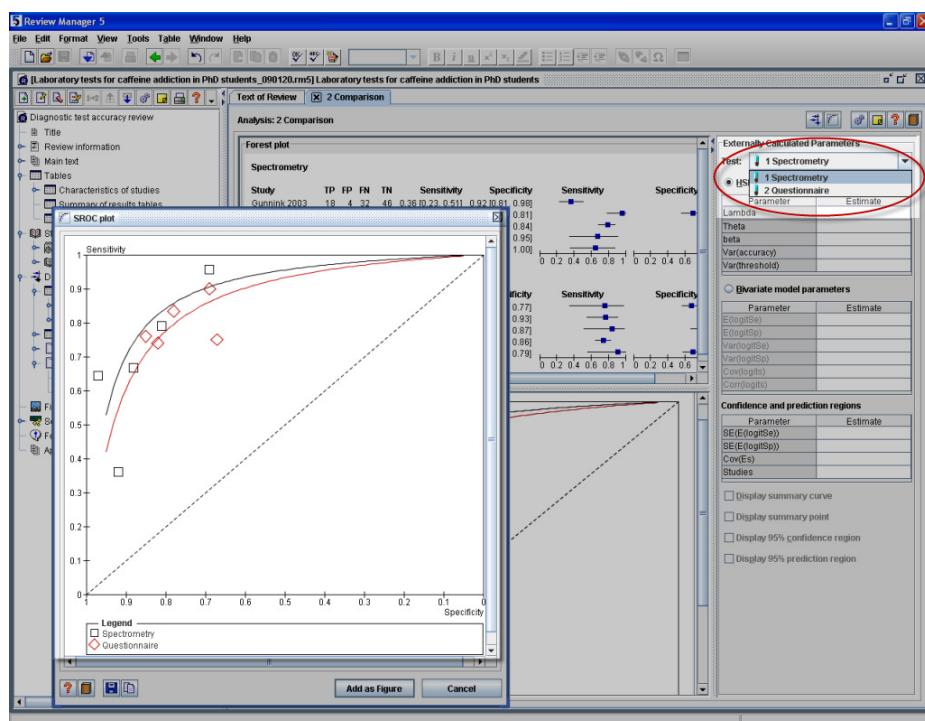

Paired Test Comparisons

RevMan enables you to restrict the Analyses to the data from studies that looked at both tests together.

1. Click the **Properties**  button in the top right corner to open the **Analysis Properties** window.
2. Select **Analyse paired data only**.
3. Click **OK**.
4. While the two forest plots display all the data for each study, the SROC plot only shows the studies that contain data for both tests, linked to each other by a solid line.
5. The curves no longer reflect the data points displayed, and should be removed.
6. Click the **Properties**  button, select the SROC plot tab and deselect **Display SROC curve(s)** option.
7. Note that in one study, the questionnaire seems to have a much higher sensitivity than spectrometry, while in the other study, spectrometry has a higher sensitivity.

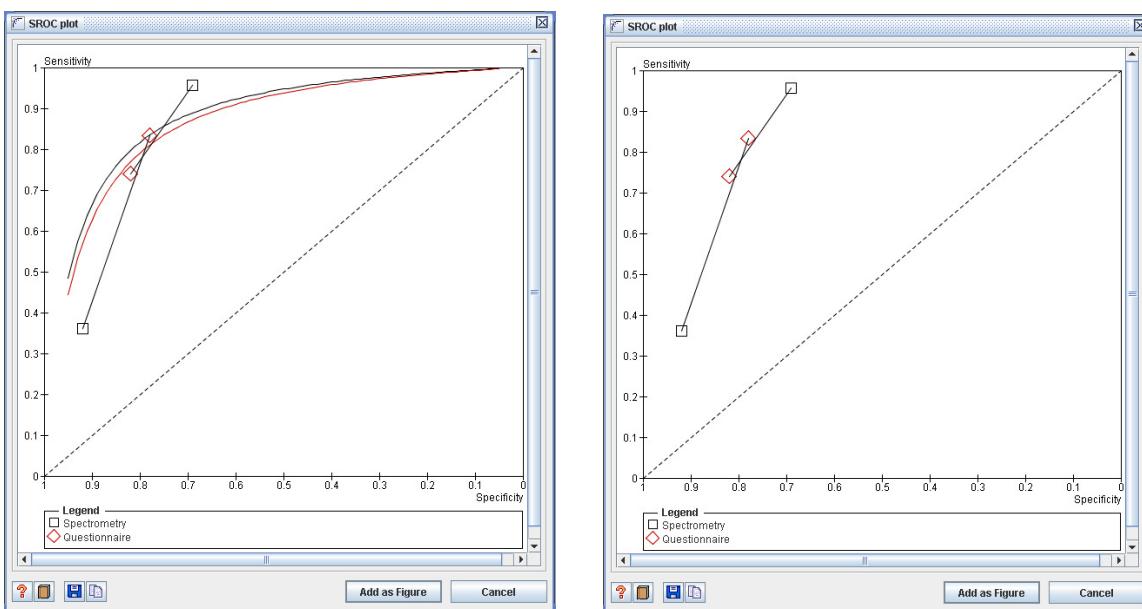


Investigating heterogeneity


RevMan allows you to investigate heterogeneity visually. You can examine whether the accuracy of studies differ by certain characteristics of studies. A stratified analysis of studies will be performed showing whether accuracy of studies differ by subgroups of studies. In RevMan no formal testing will be performed to determine whether results are significantly different. RevMan will only show a summary ROC curve for each subgroup. The forest plots may show the values of both categorical and continuous covariates, whereas the summary ROC plots will only be plotted for categorical covariates.

1. In *outline pane*, select **Analyses**.
2. Click the **Add Analysis**  button on the *outline pane* toolbar. This will open the **New Analysis Wizard**.

3. Enter the **Name** ‘Effect of questionnaire type’ and click **Next**.
4. Select **Investigate sources of heterogeneity**, then select **Questionnaires**, click **Next**, click **Finish**.

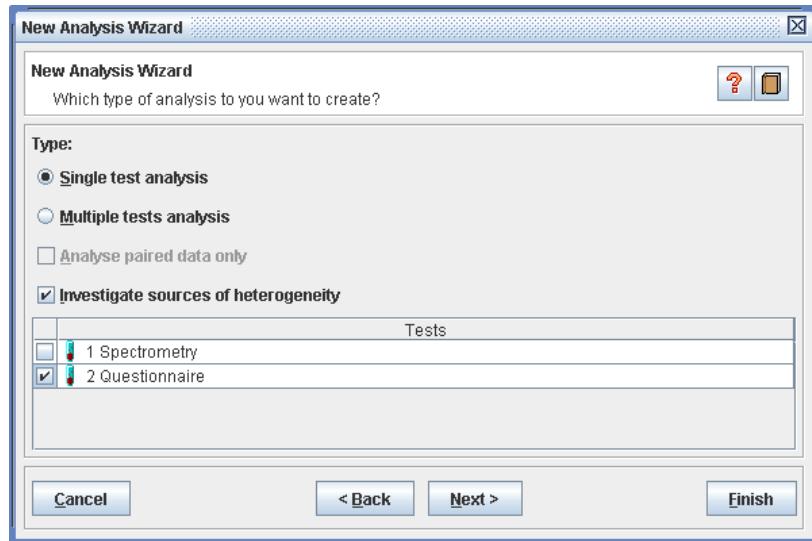

5. After clicking finish, the Analysis pane will appear, although at first there will be no indication that heterogeneity has been taken into account. To address the effect of different questionnaire types, click the **Properties** button in the top right corner. This will open the **Analysis Properties** window.
6. Open the tab **Sources of heterogeneity** and select **Covariate**.
7. Click on the drop down menu and select **Questionnaire Type**.
8. Click **OK**.
9. In the **Externally Calculated Parameters** area, a drop-down menu appears to enable you to enter the externally calculated parameters for each covariate separately.

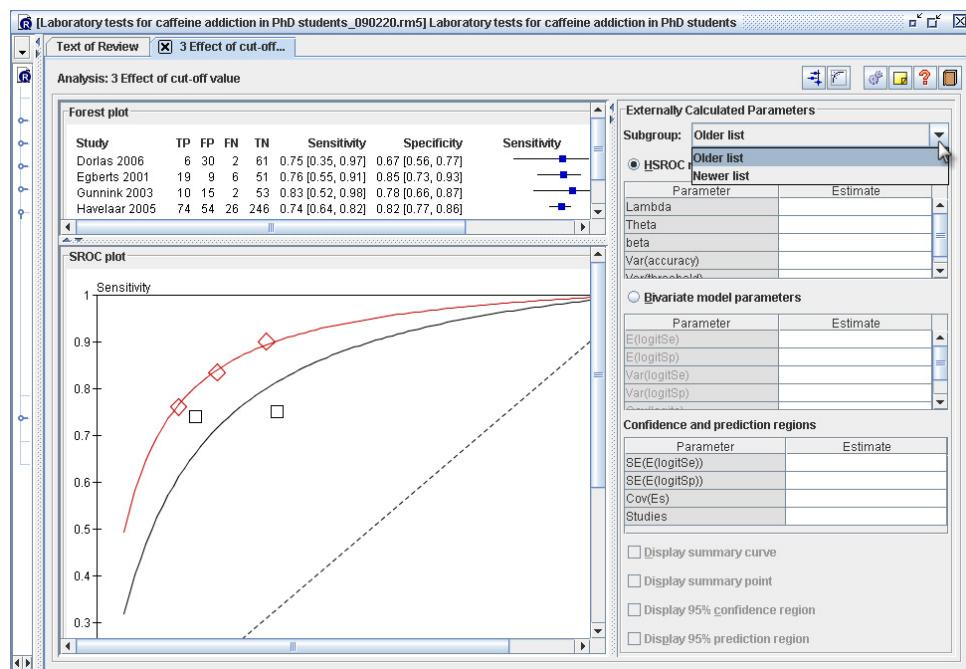

Adding a figure to your review

A small number of key summary ROC plots and forest plots can be included as figures in the text of your review, which will be published as a PDF. Forest plots for all outcomes will still be accessible as additional data in (the full version of) your published review.

1. In the *outline pane*, double-click the **Spectrometry** analysis to open the **Spectrometry** tab in the *content pane*.
2. Click the **Forest plot**  button or the SROC plot  in the top right corner of the **Spectrometry** tab. This will display the plots as they will appear in the published review.
3. Click **Add as Figure**. This will open a new tab showing the forest plot, labelled as **Figure 1**.
4. In the *outline pane*, click the key icon  next to the **Figures** heading to see that the figure has been added.
5. In the *outline pane*, click on the **Spectrometry** data table to return to the Data Tables.
6. In the **Gunnink** study, change the number of true positives from 18 to 80.
7. Click on the **Figure 1** tab again. Note that the numbers have been automatically updated.
8. Click the **Save As**  button in the top right corner of the **Figure 1** tab. This allows you to save the figure as a graphic file on your computer. Click **Cancel**.
9. Click the **Text of Review** tab. Scroll up to the **Results** section of your review.
10. Click the  icon to open the **Results** section, and then the **Findings** section.
11. Click to place the cursor in the **Findings** section.
12. Click the **Insert Link**  button on the toolbar.
13. From the **Insert Link To** list, select **Figure**.
14. Select **Figure 1** and click **OK**. A link to Figure 1 will now appear in the text of the review. In the published version, a thumbnail image will appear.

Methodological quality graphs and summaries

1. In the *outline pane*, select **Figures**.
2. Click the **Add Figure**  button on the *outline pane* toolbar. This will open the **New Figure Wizard**.
3. Select **Methodological quality graph** and click **Next**.
4. Click **Finish**. RevMan will create a graphical representation of the judgments (yes, no and unclear) entered in your **Assessment of methodological quality tables**.
5. As you will see, there are no yellow bars to indicate an 'unclear' judgment. This is because no reasons for our judgments have been entered in the tables. RevMan interprets an unclear judgment with no reason given as "not applicable", and removes these judgments from the tables in your published review.
6. Repeat steps 1 and 2. This time select **Methodological quality summary**.

7. Click **Next** and then **Finish**. RevMan will create an alternative graphical representation of the judgments in the **Methodological quality tables**. You can now add links to either of these figures in the text of the review.

PART 6 – Finishing the review

Summary of results tables

When your review is complete, you should add one or more tables to briefly summarise your results. There is no fixed format for these tables in reviews of diagnostic test accuracy. For more information on how to complete a summary of findings table, see the *Cochrane Handbook for Systematic Reviews of Diagnostic Test Accuracy*, available from <http://srtda.cochrane.org>.

1. In the *outline pane*, use the key icon  to expand the **Tables** heading, and select **Summary of results** tables.
2. Click the **Add Summary of Results Table**  button on the *outline pane* toolbar. This will open a 3 by 3 empty table, free to adjust as you like.
3. You can add rows and columns to this table by right-clicking on the table.
4. Select the top row of the table, right-click and select **Merge Cells**. A merged top row can be used to state the review question and other information needed to use the table as a stand-alone document.
5. The **Summary of results** table will be added to your review.

Adding an appendix


You can include more detailed information that is not essential to the main body of your review in an appendix, such as detailed search strategies or statistical methods.

1. In the *outline pane*, select the **Appendices** heading (the last heading on the list).
2. Click the **Add Appendix** button  on the *outline pane* toolbar. This will open the **New Appendix Wizard**.
3. Enter the **Title** ‘Detailed search strategies’.
4. Click **Finish**. In the *content pane*, the Appendices section will open and your appendix title will be shown. You can now enter text or tables in the appendix.

Validating a review

Before publication, every review should be validated to make sure all sections have been completed and meet the requirements of The Cochrane Collaboration.

1. From the **Tools** menu, select **Validate as You Type**, or click the **Validate as You Type**  button on the toolbar. In the *outline pane*, any sections with problems identified will be underlined.
2. Hover the mouse over a section underlined in red. A **Validation error** message will appear. Validation errors must be corrected before publication.
3. Hover the mouse over a section underlined in blue. A **Validation warning** message will appear. It is not compulsory to correct Validation warnings.
4. From the **File** menu, select **Reports**, and then **Validation report**. A report will be displayed listing all the **Validation errors** and **Validation warnings**. You can choose to copy, print or save this report.
5. Click **OK**.

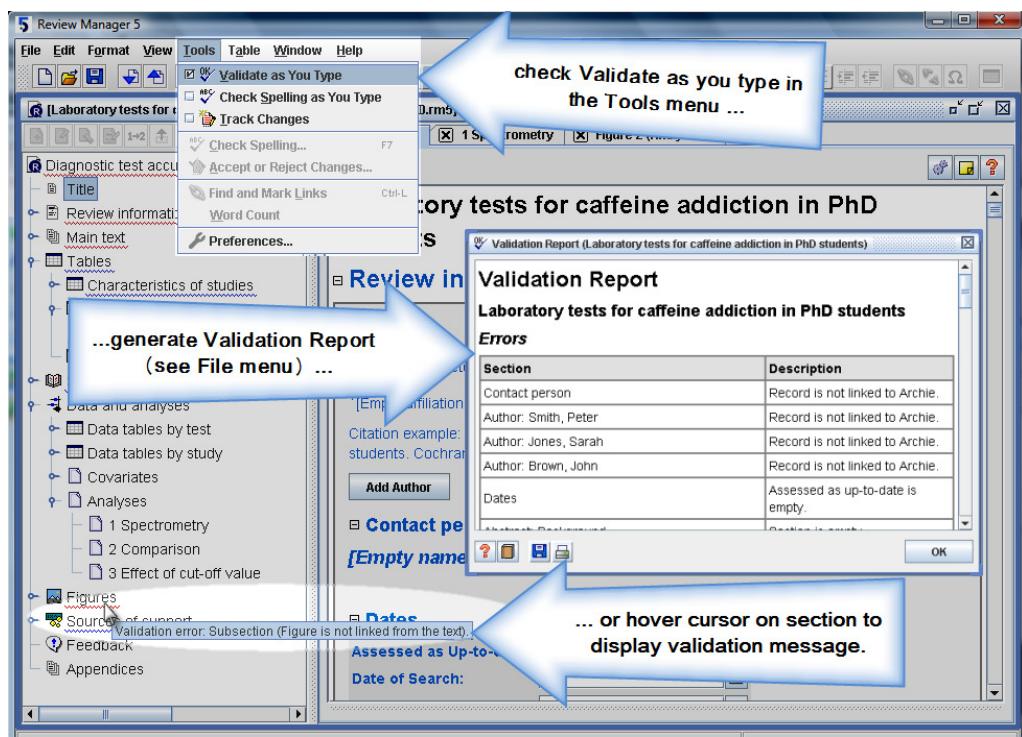

Submitting a review for editorial approval

When you have completed all the remaining sections of the review, including the **Abstract**, you will need to submit the review to your Cochrane Review Group for editorial approval. You can select this option when checking the review into the Archie database, and your Managing Editor will receive an automatic message advising that the review is ready.

You will not be able to access your review during the editorial process.

When you have made all the changes requested and the review is ready to publish, you and all your co-authors will need to sign the **Licence for Publication Form**, available from the **File** menu under **Reports**.

Note: This tutorial has been prepared by Mariska Leeflang, with assistance from the Australasian Cochrane Centre and the UK and Continental European Support Unit for diagnostic test accuracy reviews, and with funding from The Cochrane Collaboration Steering Group.